



Beautiful Access,
Practical Isolation, in a
**Comfortable
Way!**





Umbrella™

TONGUE, LIP, & CHEEK RETRACTOR

The new Umbrella retractor was developed to give you a clear treatment field while making patient comfort a top priority. See the benefits that go beyond a simple retraction of the soft tissues:



- Easy to place
- Disposable
- Naturally and gently helps the patient hold their mouth open without pulling or stretching their lips
- A new, innovative tongue-retraction design allows the tongue to comfortably rest behind the tongue guard, keeping it back and away from the working area
- Designed with anatomically placed/shaped bumpers, so clinicians can rest a hand on the patient's mouth without causing discomfort
- Provides relief and comfort to gaggers—it doesn't initiate the gag reflex for most
- Can be kept in place when checking bite
- Not sure what size your patient needs? Err on the side of Large to ensure adequate retraction.




STEP 1.

Simply press the tabs on the Umbrella tongue, lip, and cheek retractor together, ensuring the arrows on the top tab are pointing up, to prepare for insertion.



Do not place the retractor upside down.

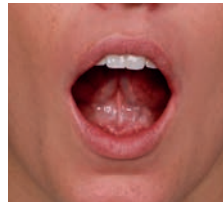


STEP 2.

Ask the patient to place the tip of their tongue on the roof of their mouth.



STEP 3.

Choose one side of the mouth in which to start, and then comfortably insert the other side of the retractor into the cheek.

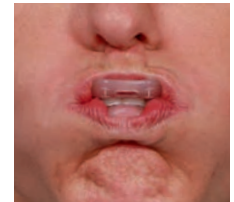

STEP 4.

Use the tabs to center the retractor with the patient's mouth.



STEP 5.

Check that the patient's tongue is resting comfortably behind the guard, ensuring easy access.



Do not place the tongue guard on top of the patient's tongue.

4870 - Umbrella Retractor Medium 5pk
5 x Umbrella Retractors

4871 - Umbrella Retractor Medium 20pk
20 x Umbrella Retractors

5162 - Umbrella Retractor Medium 40pk
40 x Umbrella Retractors

5256 - Umbrella Retractor Large 5pk
5 x Umbrella Retractors

5257 - Umbrella Retractor Large 20pk
20 x Umbrella Retractors

5258 - Umbrella Retractor Large 40pk
40 x Umbrella Retractors



SCAN
TO VIEW
PRICING



© 2024 Ultradent Products, Inc. All rights reserved.
1009272-PKAR01 102723

800.552.5512 | ULTRADENT.COM

