

MODERN CEMENTING TECHNIQUE - HIP



IMPLANT

CEMENT

BONE



ZIMMER BIOMET
Your progress. Our promise.®

Modern Cementing Technique Hip for improved clinical outcome^{1,2,3}

Modern Cementing Technique Hip is a documented and clinically proven procedure.^{1,2} The objective is to improve mechanical interlock between bone and cement in order to establish a durable interface. The use of distal plug, cement gun, pulsatile lavage and cement pressurizing devices have been shown to significantly improve long-term outcomes.³ Modern Cementing Technique, compared to earlier techniques, has been linked to a 20% reduction for the risk of revision.⁴

Zimmer Biomet offers a comprehensive portfolio of products and educational courses to supporting the key elements of Modern Cementing Technique:

Bone Bed Preparation

Preparation of the bone bed with a high pressure pulse lavage system enhance the cement penetration and the interface strength.⁵⁻⁷ It also reduces the risk of fat embolism.⁵

Bone Cement

PMMA bone cements fill the space between prostheses and bone. This helps to transmit and evenly distribute loads.

Mixing and Delivery

Mixing and collecting the cement under vacuum reduces both micro and macro pores to a minimum, thereby increasing fatigue life.^{3,8-11}

Pressurization

Pressurization improves interface between bone and cement¹⁰





IMPLANT

CEMENT

BONE



Modern Cementing Technique Hip Key Surgical Steps

Acetabulum

Bone Bed Preparation

Transversal ligament

Remove central osteophytes and be sure to save the transversal ligament in order to obtain optimal pressurization when the cup is cemented.



Reaming

Ream the acetabulum at the anatomical site to the size determined at preoperative planning. Save the part of the subcondral bone. The reaming should be over-sized in order to obtain a 2-3 mm cement mantle.




Anchorage Holes

In order to remove as little bone as possible, drill or impact five to eight holes, six mm deep, in the cranial and central parts of the acetabulum. Additional smaller diameter holes can be added in the most cranial parts. These anchorage holes increase the contact area between bone and cement, providing for better fixation.




Pulse Lavage

In order to obtain micro-interlock, use high-pressure pulse lavage and suction repeatedly. Always use a nozzle with front orifices and a splash shield. Clean repeatedly until clear fluid is received in the return line to reduce the amount of debris, blood, bacteria and fat.⁵⁻⁷

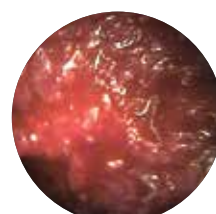

Before

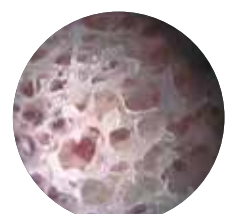

After

Vacuum Mixing of Bone Cement

The cement is mixed and collected in the cartridge under vacuum. The cartridge is then positioned in the cement gun.




Bone Cement Delivery

Break off the long nozzle at the circular mark.



Deliver the doughy bone cement into the acetabulum using the cement gun with short nozzle.

In order to minimize lamination of the blood, apply the vacuum mixed bone cement as a bolus. This requires high viscosity bone cement. Delivery time varies depending on type of cement used.

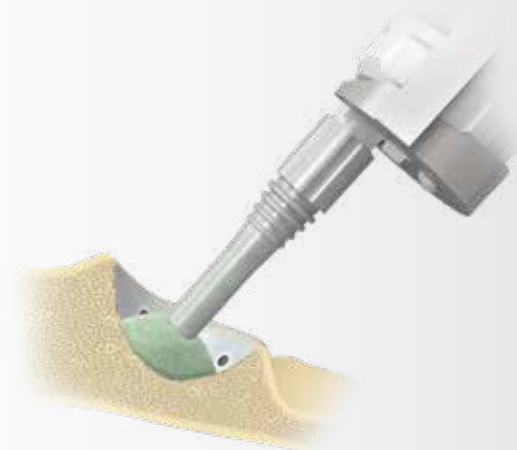



ⓘ **Note:** Protect your hands and prevent contact allergy by using an extra pair of PE gloves. Latex gloves provide no protection from monomer penetration.

Pressurization

Pressurization

Pressure is applied immediately by using the acetabular pressurizer. Maintain pressure until the bone cement is sufficiently doughy.


Collected data has shown that high pressurization is needed to achieve micro-interlock.² This reaffirms the necessity of working in a contained cavity and not removing the transversal ligament.

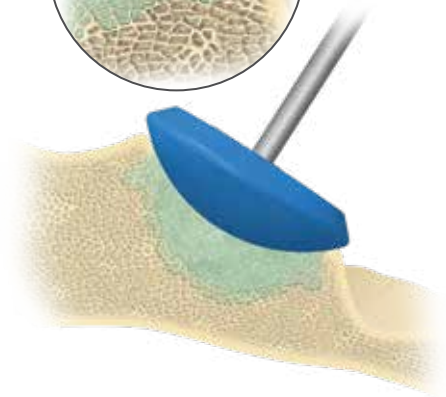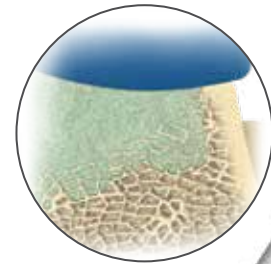

Even Cement Mantle


You should have a 2-3 mm cement mantle around the cup.

An even cement mantle means better stress distribution and reduces risk of cement mantle failure.

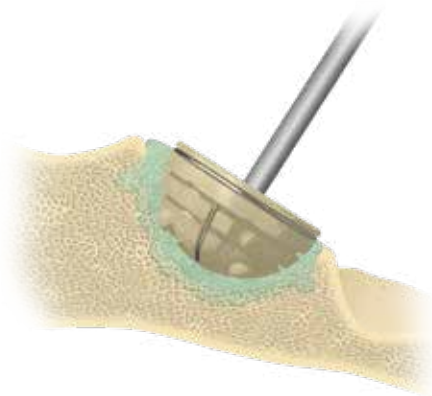

Protect the Cup

Use the cup protector in order to prevent cement contamination in the articulating surface.

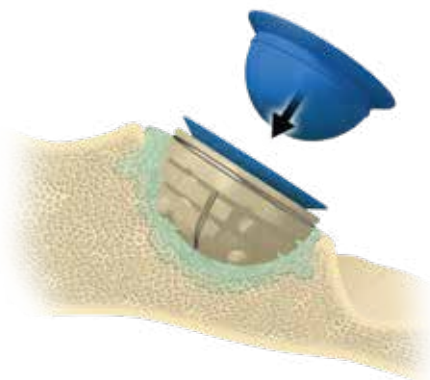

Femur

Bone Bed Preparation

Reaming

Use a straight reamer to open up the femoral canal. Be sure to open both the posterior and the lateral walls sufficiently.

Insufficient reaming will make it impossible to align the femoral stem properly in both planes. Start with the smallest available rasp impactor and increase until the size determined at preoperative planning has been achieved. If necessary, be prepared to use a flexible reamer to enlarge the isthmus.



Distal Plug

Determine the plug size by using the special instrument or a flexible reamer. Insert the appropriate distal femoral plug. Place the plug 10 mm distal to the intended level of the tip of the stem. Resorbable biodegradable plugs make the removal during revision surgery unnecessary.



Pulse Lavage

Use a nozzle with side orifices, allowing the pulse lavage to act perpendicular to the bone surface.

Clean repeatedly until clear fluid is received in the return line to reduce the amount of debris, blood, bacteria and fat.⁵⁻⁷



Vacuum Mixing of Bone Cement

The bone cement is mixed and collected in the cartridge under vacuum.

The cartridge is then positioned in the cement gun.



Pulse Lavage

Make a final pulse lavage before injecting the cement.

Clean repeatedly until clear fluid is received in the return line to reduce the amount of debris, blood, bacteria and fat.⁵⁻⁷



Bone Cement Delivery

The doughy bone cement is delivered in retrograde fashion.




Pressurization

Break off the long nozzle at the circular mark.




Place the support plate and the femoral pressurizer on the short nozzle.

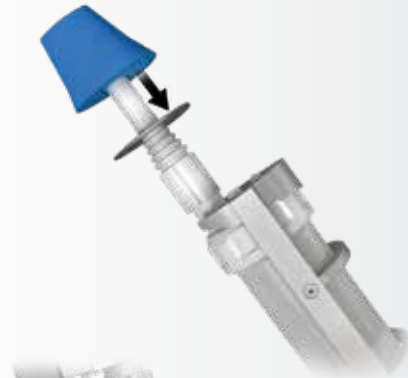


Apply the proximal seal and pressurize the bone cement. A positive sign of pressurization is marrow extrusion in the greater trochanter (the so-called sweating trochanter sign).

As with the acetabulum, maintain pressure until the bone cement is sufficiently doughy to withstand bleeding. The time varies depending on the type of cement used.

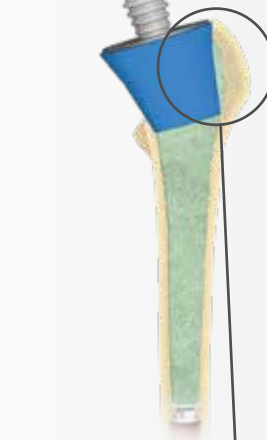

Bone cement penetration into cancellous bone insert.

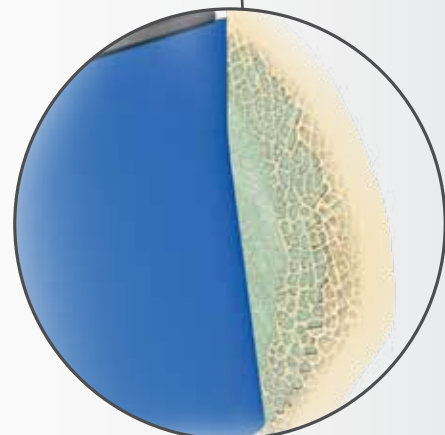

Pressurization

Distal centralizer

Apply the distal centralizer to the stem. The centralizer should be used to avoid varus or valgus malposition of the stem in both planes.



Introduction of the stem

Gently introduce the stem and hold it in place until the cement has polymerized.

Even cement mantle

With the stem in final position, you should have a 2–3 mm cement mantle around the stem and approximately 10 mm between the tip of the stem and the plug. This will yield optimal stress distribution. Finally, reduce the hip.



Modern Cementing Technique Hip

The method for improved clinical outcome.^{1,2,3}

The modern total hip replacement has effectively helped to improve quality of life for patients with osteoarthritis of the hip.^{2,2} Modern and reproducible cementing technique has definite advantages for our patients as well as for the health care providers.³

Long Term Implant Survival

Good implant stability

Micro-interlock


Modern Cementing Technique Hip

Good initial stability of the implant is one important factor to obtain good and long lasting results. Each step of Modern Cementing has been associated with approximately 20% reduction of revision for aseptic loosening.⁷

Micro-interlock and a homogenous bone cement mantel around the implant are important factors in achieving good implant stability

Micro-interlock is the mechanical interlock between bone and cement. It is important to have an optimal micro-interlock in order to achieve a durable interface.

The goal of Modern Cementing Technique in general is to improve both the long-term properties of the bone cement itself and the cement-bone interface.

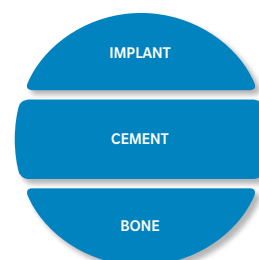

National Hip Registers

Since 1979, the Swedish Total Hip Replacement Registry has documented cement techniques. More than 400 000 surgeries have been recorded and the data demonstrates significantly improved longevity of hip implants when Modern Cement Technique is utilized.² The significant variables are distal plug, pulsatile lavage, proximal seal and vacuum mixing of the bone cement. Each of these steps has been associated with approximately 20% reduction of revision for aseptic loosening.⁷

Cement-Bone Interface

Bone Bed Preparation

Preparation of the bone bed with a pulsative lavage system, like the Pulsavac® Plus Wound Debridement System, helps to ensure solid cement fixation.¹

- To obtain proper cement penetration and fixation into the cancellous bone¹⁴
- Reduce the risk for revision due to aseptic loosening¹
- Reduce the risk for fat embolism and significantly minimize circulatory changes.^{5,16}



Bone Cement


Polymethyl methacrylate (PMMA) bone cements fill the space between prostheses and bone, transmitting and evenly distributing loads. Antibiotic-loaded bone cement reduces the risk of infection in both primary and revision operations.^{2,13}

High viscosity cements have been shown to offer a lower incidence of revision and aseptic loosening in total hip replacement.^{2,13}

The handling properties of the bone cement are highly dependent on the temperatures of the cement and the operating room. Higher temperatures make for a shorter working phase and a faster setting time.

High viscosity bone cement like Refobacin® Bone Cement R can be pre-chilled if a longer working phase is required. The main considerations are:

- Good mechanical properties²¹
- Consistent handling properties²²
- High antibiotic release

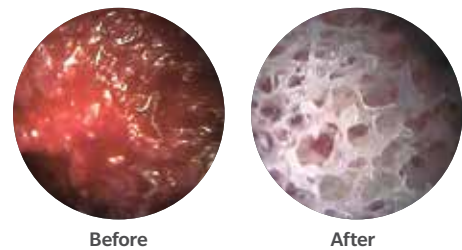


Before

After




Factors Influencing Bone Cement Handling Characteristics



Temperature

- Cement temperature
- Body temperature
- OR temperature
- Storage temperature


Higher temperature = faster setting time



Time

- From storage to mixing
- From mixing to application
- Application time
- Setting time
- Start the clock immediately when the powder and liquid meet


Test – get to know your cement



Mixing Method

- Vacuum (closed system)
- Bowl (open air)


Read relevant information



Type of Cement

- Low viscosity
- Medium viscosity
- High viscosity


Read relevant information

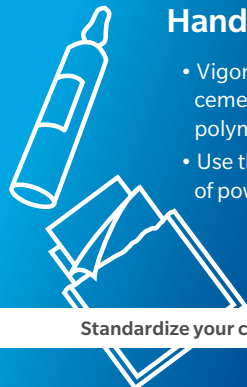

Handling

- Vigorous kneading of the cement may speed up the polymerization process
- Use the right proportion of powder and liquid

Standardize your cement handling



Improved Working Environment

- Use an extra pair of PE gloves
- A closed system minimizes monomer fumes and skin contact
- Use a combined system for mixing and delivery

Rubber gloves do not protect against monomer

Mixing and Collection Under Vacuum

Mixing and collecting the bone cement under vacuum reduces both micro and macro pores.^{3,8-11}

- Improved cement strength and fatigue life^{3,17}
- Lower risk of aseptic loosening due to cracks^{9-11,17}
- Delivery of reproducible results^{21,22}
- Less exposure to monomer fumes^{18,19}

Delivery


- Acetabulum: Deliver the doughy bone cement into the acetabulum using the cement gun with the short nozzle. It is recommended to apply the cement as a bolus.
- Femur: Deliver the doughy bone cement into the femur using the cement gun with the long nozzle. The doughy cement should be delivered in a retrograde fashion. The Swedish Hip Register has shown that retrograde cement filling in the femur reduces the risk of revision.²

Pressurization


High pressure increases penetration into cancellous bone to achieve micro-interlock

- Reduces cement porosity and increase fatigue strength^{3,12}
- Improved interface between bone and cement^{3,14}
- Better stress distribution³

Cement mixed under vacuum^{9*}

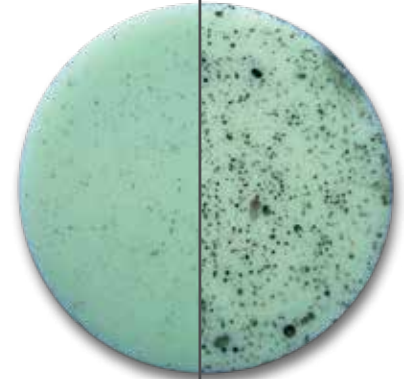

Cement mixed at atmospheric pressure^{9*}

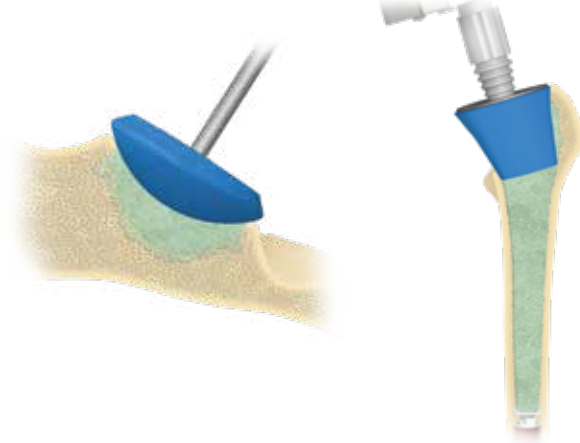

* Lab test results not necessarily indicative of clinical performance.

Solutions for Modern Cementing Technique Hip

Optipac® Vacuum Mixing System

A closed vacuum mixing system, pre-packed with bone cement

Proven*, Strong, Safety

- On the market since 2008 and based on Optivac® Technology since 1993
- Designed to improve cement fatigue life by mixing and collection under vacuum^{1,8-10}
- Featuring SoftPac™ Technology ensures no breaking of glass ampoules.
- Standardized and reproducible method, improved bone cement quality⁹
- Minimizing exposure to monomer fumes^{18,19}
- No direct contact with bone cement during mixing and delivery
- No breaking of glass ampoules, no risk of glass injury



Optipac | SINCE 2008

Optivac® Vacuum Mixing System

Designed for mixing and collection under vacuum,

Optivac Vacuum Mixing System reduces both microporosity and macroporosity⁸⁻¹¹

- Designed to improve cement strength and fatigue life⁸⁻¹¹
- On the market since 1993
- Unmatched in documentation^{8-11, 18 **}



Refobacin® Bone Cement R

- High viscosity cement
- Green color of cement allows for easy recognition during surgery
- Reliable mechanical performance based on international standard laboratory testing^{20,21*}




*Technology tried in practice of mixing and collecting bone cement in a system under vacuum.


**Laboratory testing is not necessarily indicative of clinical performance. In published articles and lab reports

Ordering Information


Optipac Refobacin® Bone Cement R

Product	Description	Part Number	Units/box
	Optipac Knee Refobacin® Bone Cement R	4709500392-3	8
	Optipac 40 Refobacin® Bone Cement R	4710500394-3	8
	Optipac 60 Refobacin® Bone Cement R	4711500396-3	8
	Optipac 80 Refobacin® Bone Cement R	4712500398-3	8
	Optipac Hip Set (40 + 80) Refobacin® Bone Cement R	4740500394-3	4


Optipac Refobacin® Plus Bone Cement

Product	Description	Part Number	Units/box
	Optipac Knee Refobacin® Plus Bone Cement	4719502082-3	8
	Optipac 40 Refobacin® Plus Bone Cement	4720502083-3	8
	Optipac 60 Refobacin® Plus Bone Cement	4721502084-3	8
	Optipac 80 Refobacin® Plus Bone Cement	4722502117-3	8


Optipac Refobacin® Revision

Product	Description	Part Number	Units/box
	Optipac 40 Refobacin® Revision	4730501163-3	8
	Optipac 80 Refobacin® Revision	4732501165-3	8


Refobacin® Bone Cement R

Product	Description	Part Number	Units/box
	Refobacin® Bone Cement R 1x20	3003920001-3	20
	Refobacin® Bone Cement R 2x20	3003920002-3	20
	Refobacin® Bone Cement R 1x40	3003940001-3	20
	Refobacin® Bone Cement R 2x40	3003940002-3	20
	Refobacin® Bone Cement R 1x60	3003960001-3	12







Refobacin® Plus Bone Cement

Product	Description	Part Number	Units/box
	Refobacin® Plus Bone Cement 1x20	3020820401-3	20
	Refobacin® Plus Bone Cement 2x20	3021180001-3	20
	Refobacin® Plus Bone Cement 1x40	3020830401-3	20
	Refobacin® Plus Bone Cement 2x40	3021170001-3	20
	Refobacin® Plus Bone Cement 1x60	3020840401-3	12

Refobacin® Revision

Product	Description	Part Number	Units/box
	Refobacin® Revision 1x40	3011630001-3	20

Mixing & Delivery and Bone Bed Preparation


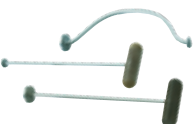
Product	Description	Part Number	Units/box
	Optivac S	4161	10
	Optivac M	4160	10
	Optivac L	4152	10
	Optivac Hip Set	4150	10
	Optivac+ Hip Set	4250	10
	Optigun™	4193	1
	Optigun Ratchet	4195	1
	Pulsavac® Plus Hip Kit	00-5150-482-00	10
	Optivac+ Hip Set	4250	10

Pressurizers

Femur Pressurizers

Product	Description	Part Number	Units/box
	Pressurizer Femur II	430900	5
	Cement Press Support Plate II	4197	1





Acetabular Pressurizers

Product	Description	Part Number	Units/box
	Acetabular Pressurizer 50 mm	4316	1
	Acetabular Pressurizer 57 mm	4317	1
	Acetabular Pressurizer 63 mm	4321	5
	Acetabular Pressurizer 71 mm	4322	1
	Handle For Acetabular Pressurizer 4321, 4322	4327	1
	Handle For Acetabular Pressurizer 4316, 4317	4318	1
	Angle Handle Acetabular Press, compatible with all acetabular pressurizers	4319	1

Ordering Information




Bone Bed Preparation

Quik-Use® Femoral Bone Prep Kits


Product	Description	Part Number	Units/box
	Quik-Use Femoral Bone Cement Prep Kit	00-5049-055-00	1
	Quik-Use Femoral Bone Cement Prep Kits	00-5049-055-10	10
	Femoral Bone Brushes	00-5059-013-00	10
	Ortho-Dri® Absorbent Packing	00-5059-080-00	5
	Quik-Use Curette	00-5049-053-00	10

Polyethylene Bone Plugs

Allen Bone Plugs

Product	Description	Part Number	Units/box
	Allen Plug HDPE/BASO4 7C/12FL	00-8011-020-12	1
	Allen Plug HDPE/BASO4 8C/16FL	00-8011-020-16	1
	Allen Plug HDPE/BASO4 10C/20FL	00-8011-020-20	1
	Allen Plug HDPE/BASO4 12C/24FL	00-8011-020-24	1
	Allen Plug HDPE/BASO4 14C/28FL	00-8011-020-28	1
	Allen Plug HDPE/BASO4,16C/32FL	00-8011-020-32	1
	12 mm/16 mm Plugs with Disposable Inserter	00-8011-010-01	1
	20 mm/24 mm Plugs with Disposable Inserter	00-8011-020-01	1
	28 mm/32 mm Plugs with Disposable Inserter	00-8011-030-01	1
	Allen Plug Inserter Metal	00-8011-210-00	1

Stühmer/Weber Bone Plugs

Product	Description	Part Number	Units/box
	Stühmer/Weber Intramed Plug W/out Drain, Size 1	0961	1
	Stühmer/Weber Intramed Plug w/out Drain, Size 2	3202	1
	Stühmer/Weber Intramed Plug w/out Drain, Size 2.5	3208	1
	Stühmer/Weber Intramed Plug w/out Drain, Size 3	3203	1
	Stühmer/Weber Intramed Plug w/out Drain, G Size 3.5	3209	1
	Stühmer/Weber Intramed Plug w/out Drain, Size 4	3204	1
	Stühmer/Weber Intramed Plug w/out Drain, Size 5	3205	1
	Stühmer/Weber Intramed Plug w/out Drain, Size 6	3206	1
	Stühmer/Weber Intramed Plug w/out Drain, Size 7	3207	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 1	0981	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 2	3222	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 2.5	3228	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 3	3223	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 3.5	3229	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 4	3224	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 5	3225	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 6	3226	1
	Stühmer/Weber Intramed Plug w/ Drain, Size 7	3227	1

Instrumentation

Product	Description	Part Number	Units/box
	Plug Inserter	5950	1
	Ø 7.1 mm Measuring Cone, Size 1.0	5951	1
	Ø 9.1 mm Measuring Cone, Size 2.0	5952	1
	Ø 9.8 mm Measuring Cone, Size 2.5	5958	1
	Ø 11.5 mm Measuring Cone, Size 3.0	5953	1
	Ø 12.7 mm Measuring Cone, Size 3.5	5959	1
	Ø 13.9 mm Measuring Cone, Size 4.0	5954	1
	Ø 16.4 mm Measuring Cone, Size 5.0	5955	1
	Ø 19 mm Measuring Cone, Size 6W	5961	1
	Ø 22 mm Measuring Cone, Size 7W	5962	1

References

1. Malchau H, et al. Prognosis of Total Hip Replacement. The National Hip Arthroplasty Register: 9-11,1996.
2. Malchau H, et al. The Swedish Total Hip Replacement Register. *The Journal of Bone and Joint Surgery*. 84A: 2-20, 2002.
3. Breusch S, et al. Optimal Cementing Technique - The Evidence: What is Modern Cementing Technique? In: The Well-Cemented Total Hip Arthroplasty. Chapter 6.1. Springer Verlag 2005.
4. Swedish Hip Arthroplasty Register. Annual Report 1998.
5. Breusch S, et al. Pulmonary Embolism in Cemented Total Hip Arthroplasty. The Well-Cemented Total Hip Arthroplasty. Heidelberg, pp. 320-331, 2005. DOI 10.1007/3-540-28924-0_43.
- 6* Anglen JO, et al. The efficacy of various irrigation solutions in removing slime-producing Staphylococcus. *Journal of Orthopaedic Trauma*. 8(5):390-6, 1994.
7. Helmers, S, et al. Efficacy of Irrigation for Removal of Particulate Debris after Cemented Total Knee Arthroplasty. *Journal of Arthroplasty*. 14(5):549-52.1999.
- 8* Wang J-S, et al. Porosity of Bone Cement Reduced by Mixing and Collecting Under Vacuum. *Acta Orthopaedica Scandinavica*. 64(2):143-6, 1993.
- 9* Wang J-S, et al. Bone Cement Porosity in Vacuum Mixing Systems, in Bone Cements and Cementing Technique, Walenkamp, Murray (Eds). Springer Verlag. 2001.
- 10* Dunne N-J, et al. Influence of the Mixing Techniques on the Physical Properties of Acrylic Bone Cement. *Biomaterials*. 22(13) 1819-26, 2001.
- 11* Wilkinson JM, et al. Effect of mixing technique on the Properties of Acrylic Bone-Cement. *The Journal of Arthroplasty*. 15:663-667, 2000.
- 12* Reading, AD, et al. A Comparison of 2 Modern Femoral Cementing Techniques: Analysis by Cement-bone Interface Pressure Measurements, Computerized Image Analysis, and Static Mechanical Testing. *The Journal of Arthroplasty*. 15(4):479-87, 2000.
13. Havelin L-I, et al. Prospective Studies of Hip Prostheses and Cements. A presentation of the Norwegian Arthroplasty Register 1987-1999. Scientific Expedition at the 67th Annual Meeting of the American Academy of Orthopaedic Surgeons. March 15-19, 2000 Orlando, USA.
- 14* Kalteis T, et al. An experimental comparison of different devices for pulsatile high-pressure lavage and their relevance to cement intrusion into cancellous bone. *Archives of Orthopaedic and Trauma Surgery*. 127(10):873-7.2007.
15. Christie J, et al. Medullary Lavage Reduces Embolic Phenomena and Cardiopulmonary Changes During Cemented Hemiarthroplasty. *The Journal of Bone and Joint Surgery*. (Br). 77-B:456-9, 1995.
16. Breusch S. Cementing techniques in Total Hip Replacement: factors influencing survival of femoral components. In: Bone Cements and Cementing Technique; Walenkamp GHIM, Murray DW (eds) Berlin, Heidelberg, Springer Verlag 2001.
- 17* Report from SP Technical Research Institute of Sweden Airborne Methyl Methacrylate Monomer During the use of Different Bone Cement Mixing Systems. (2007 08 13).
- 18* Schlegel UJ, et al. Pre-Packed Vacuum Bone Cement Mixing Systems. A Further Step in Reducing Methylmethacrylate Exposure in Surgery. *Annals of Occupational Hygiene*.40(8) 955-61, 2010.
19. ISO 5833, Implants for Surgery – Acrylic Resin Cements (2002).
- 20* Data on file at Zimmer Biomet, internal laboratory testing results, Aug 2016. Laboratory testing is not necessarily indicative of clinical performance.
- 21* Data on file at Zimmer Biomet. Internal laboratory testing results of Fatigue performance. Project Number: BC0051 - Part A&D; Optipac Refobacin Bone Cement R. Sept 2017. Laboratory testing is not necessarily indicative of clinical performance.
22. Learmonth, D.I. et al. The Operation of the Century: Total Hip Replacement. *The Lancet*. 370:1508-19, 2007.

*Lab test results not necessarily indicative of clinical performance.

All content herein is protected by copyright, trademarks and other intellectual property rights, as applicable, owned by or licensed to Zimmer Biomet or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet.

This material is intended for health care professionals. Distribution to any other recipient is prohibited.

For product information, including indications, contraindications, warnings, precautions, potential adverse effects and patient counseling information, see the package insert and www.zimmerbiomet.com.

Check for country product clearances and reference product specific instructions for use. Not for distribution in France.

Zimmer Biomet does not practice medicine. Each surgeon should exercise his or her own independent judgment in the diagnosis and treatment of an individual patient, and this information does not purport to replace the comprehensive training surgeons have received. As with all surgical procedures, the technique used in each case will depend on the surgeon's medical judgment as the best treatment for each patient. Results will vary based on health, weight, activity and other variables. Not all patients are candidates for this product and/or procedure. Caution: Federal (USA) law restricts this device to sale by or on the order of a surgeon. Rx only.

Refobacin® is a trademark of Merck KGaA.

© 2019 Zimmer Biomet



 **Legal Manufacturer**
Biomet France
Plateau de Lautagne
26000 Valence
France
Phone: 0033 4 75 75 91 00