# BRAUN

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21332986/VI-25

[www.braun.com](http://www.braun.com)

PL1

## Braun Infolines

English	4	<b>UK</b>	<b>0800 783 7010</b>
Polski	11	<b>IE</b>	<b>1 800 509 448</b>
Český	19	<b>PL</b>	<b>801 127 286</b>
Slovenský	26	<b>CZ</b>	<b>800 113 322</b>
Magyar	33	<b>SK</b>	<b>800 333 233</b>
Hrvatski	40	<b>HU</b>	<b>(06-1) 451-1256</b>
Slovenščina	47	<b>HR</b>	<b>01/ 6690 330</b>
Türkçe	54	<b>SI</b>	<b>802 822</b>
Română	62	<b>TR</b>	<b>0850 220 0911</b>
Ελληνικά	70	<b>RO</b>	<b>0800 890 194</b>
Български	78	<b>GR</b>	<b>800 801 3458</b>
Українська	86	<b>BG</b>	<b>008001184973</b>
		<b>UA</b>	<b>0 800 505 000</b>


Type 6034

Braun GmbH  
Sulzbacher Straße 40  
65824 Schwalbach / Germany  
[www.braun.com](http://www.braun.com)  
[www.service.braun.com](http://www.service.braun.com)

CyDen Ltd., Block A,  
Bay Studios Business Park,  
Fabian Way, Swansea SA18QB, UK
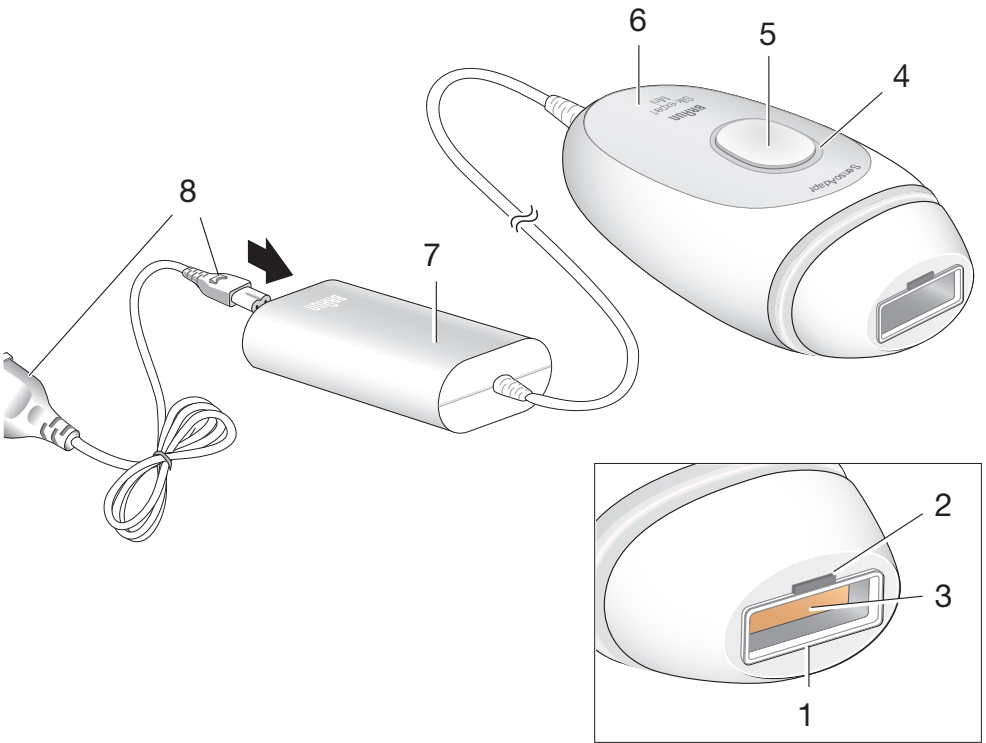
## Skin tones and body hair types

✓ = suitable / odpowiedni / vhodné / vhodné / megfelelő / prikladno / Kullanim için uygun / primerno / adecvat / κατάλληλο / подходящий / підходить

✗ = not suitable / nieodpowiedni / použití není vhodné / nevhodné / nem megfelelő / nije prikladno / kullanim için uygun değil / ni primerno / neadecvat / ακατάλληλο / не e подходящий / не підходить

# English

Our products are designed to meet the highest standards of quality, functionality and design. Before using the device, please read the use instructions carefully, including all contraindications, warnings and safety information, and keep them for future reference.

Using professional light technology, Braun Silk-expert Mini Intense Pulsed Light (IPL) enables long-term hair growth reduction after a certain period of treatment.

Ideal body areas include:

Women	Men
<p>Legs, arms, underarms, intimate area*, face (below cheek bones), fingers, toes</p> <p><b>CAUTION DO NOT USE...</b> around the eyes, along the forehead, on scalp, nipples</p>	<p>Shoulders and below (chest, back, underarms, arms, stomach, intimate area*, legs, fingers, toes)</p> <p><b>CAUTION DO NOT USE...</b> on face, neck, scalp, nipples</p>
<p><b>*Intimate area: VIO</b></p> <div style="display: flex; justify-content: space-around; align-items: center;"> <div style="text-align: center;">  <p>Bikini area</p> </div> <div style="text-align: center;">  <p>between legs</p> </div> <div style="text-align: center;">  <p>buttocks and skin around anus (peri-anal)</p> </div> </div> <p><b>DO NOT USE...</b> in mucous membranes such as labia minora, vagina or on the anus</p>	<p><b>*Intimate area:</b></p> <div style="display: flex; justify-content: space-around; align-items: center;"> <div style="text-align: center;">  <p>Area above penis</p> </div> <div style="text-align: center;">  <p>buttocks and skin around anus (peri-anal)</p> </div> </div> <p><b>DO NOT USE...</b> on penile shaft, scrotum or anus</p> <div style="text-align: right; margin-top: 10px;">  </div>

IPL is not suitable for everyone. It is most effective on light to medium skin tones and dark hair colors. Please, refer to the skin tone and hair color chart on page 3 to see if the device is suitable for you.

## Description

- 1 Treatment window
- 2 SensoAdapt™ skin tone sensor
- 3 Glass filter
- 4 Power indicator light
- 5 Treatment button
- 6 Handle
- 7 Base station
- 8 Mains cord

## How does it work?

Braun Silk-expert Mini works beneath the skin surface by targeting the melanin in the hair follicle, helping to break the cycle of hair regrowth. With continued, periodic use, the light energy helps prevent unwanted hair from reaching the skin surface.

The intelligent SensoAdapt™ skin sensor (2) continuously reads your skin tone before every flash, automatically adapting the light intensity for best efficacy and safety. This allows treatment of appropriate skin tones, but prevents treatment if your skin tone is too dark.

## Contraindications

### Do not use the device ...

- if your skin is darker than skin tone 5 shown on the skin tone chart on page 3. Since dark skin absorbs more light energy, treating very dark skin may cause discomfort/pain and adverse effects (e.g., burns, blisters, discoloration or scarring) and injure your skin.
- if you are pregnant, lactating, or under the age of 18 as this device has not been tested with these individuals.
- if you have a history of skin cancer or precancerous lesions (e.g., nevi or a large number of moles).
- if your skin has been artificially or naturally tanned recently. Your skin may be sensitive following sun exposure and particularly susceptible to side effects of IPL treatment (e.g. burn, blister, discolour or scar to your skin). Also, avoid unprotected exposure to direct sunlight after your IPL treatment.

### Certain conditions can also prohibit the use of this device.

#### Do not use the device if any of the following applies:

- You have a chronic medical condition.
- You have a chronic skin disease (e.g. psoriasis or vitiligo) or damaged skin (e.g. sunburn, cuts, open wounds or active infections) in the areas you wish to treat.
- You have varicose veins in the area you wish to treat.
- You have a known sensitivity to sunlight (photosensitivity) or are taking medication that may make the skin more sensitive (e.g. Retin A, Accutane and/or other topical retinoids).
- You have had a professional skin peel treatment on the area you wish to treat.
- You have any of the above, this device may injure your skin or make existing conditions worse. You may experience side effects such as burns, blisters and skin color changes (hyper- or hypo-pigmentation) or scarring.

If you are unsure whether this device is safe for you to use, please consult your physician or dermatologist. If you take any kind of medication regularly and over a long period of time, please consult your physician for any potential impact on skin sensitivity to light.

## Warnings

### Do not use on ...

- scalp, eyebrows or anywhere near or around the eyes. Use in these areas may permanently damage your eyes.
- a tattoo or permanent makeup in the treatment area.
- nipples, genitals or around the anus. These areas may have a darker skin color and/or greater hair density and using the device in these areas may cause discomfort/pain or injure your skin.
- dark brown or black spots such as freckles, birthmarks, moles or warts in the treatment area.
- beard / facial hair for men.

Usage on male beards or facial hair may lead to permanent or uneven results which may produce changes in facial characteristics that may not be desirable.

It is not recommended to use this device on areas where dermal fillers are present.

**Do not use** more than once a week on the same area. Damage of skin may occur after prolonged or repeated treatment on one site.

## Eye safety

**Do not look directly into the glass filter** (3) or attempt to activate the device toward the eyes. To protect the safety of your eyes and those of individuals nearby, the treatment window (1) is provided with a safety system:

The device will only flash when the entire front surface is placed against your skin.

If at any time the device flashes when the front surface is NOT in contact with the skin, STOP use immediately and contact your retailer.

Possible eye injury (potentially leading to loss of vision) or skin injury if instructions are not followed. Protect the eye from exposure. Do not use the appliance over the eye lids or close to the eye.

## Precautions

- The device is intended for a single user.
- It is not intended for use by people under 18 years.
- People with reduced physical, sensory or mental capabilities or lack of experience and knowledge, may use the device if they have been given supervision or instruction concerning the safe use of the device and understand the hazards involved.
- Children should be supervised to ensure that they do not play with the device. Always keep the device out of children's reach. Keep the device out of reach of pets and pests which could damage the device.

## Electrical and fire safety

Regularly check the device and cords for visible signs of damage. In case of damage or cracks, stop using the device. Electric shock may result from using a damaged device. In case of defects, have the device checked and repaired by an authorized service center. Unqualified repair work can lead to extreme hazards for the user. Always unplug the device after use.

### DANGER

Never attempt to repair and open the device, as this may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury.

**Do not** use in or near water and wet surfaces such as bathtubs, showers, basins or other vessels that may contain water. If the device becomes wet, it should not be used. Electric shock may result from using a wet device.

**Do not** use if the device is damaged (e.g. cracked housing, cord damage etc). Internal parts contain electric charges that are hazardous.

Regularly inspect the treatment window (1) for damage:

**Do not** use if the glass filter (3) is cracked or missing.

**Do not** use if the built-in reflector is missing or not in correct place.

**Do not** touch the glass filter during or directly after use since it will get hot.

**Do not** use petroleum-based or flammable cleaning agents (alcohol, acetones) on the device as these pose a fire risk.

**Do not** use abrasive cleaning agents. These may damage the device or scratch the glass filter.

**Do not** attempt to use this device other than for removing unwanted hair.

**Do not** override the safety mechanism inherent to the device.

**Do not** use the device on a surface other than skin. Do not press the treatment button when the treatment window (1) is not placed in full contact with the skin.

**Do not** block the cooling vents. The device may overheat and pose an electric shock risk or fire hazard.

## How to use

### Preparing the treatment area

Remove any visible hair on the skin before using the device. You may use your preferred hair removal method. To keep the treatment window clean and intact and prevent hair burns, ensure there is no hair remaining on the skin surface. Clean the treatment area and pat it dry.

### To power on the device:

Use the mains cord (8) to connect the base station to an electrical outlet. Press the treatment button (5) to activate the device. The power indicator light (4) lights up once and the fan starts running. This means the device is in READY mode. During operation, the device will heat up slightly.

### Automatic light adaptation

When the skin tone sensor (2) has skin contact, the device reads your skin tone and selects the appropriate light intensity which is displayed by the power indicator light color.

Pink:	high intensity for light skin tones
Light pink:	medium intensity for medium skin tones

White:	low intensity for dark skin tones
Red:	Invalid skin tone / device will not flash

After 10 minutes without use, the device is automatically deactivated. To reactivate it, press the treatment button.

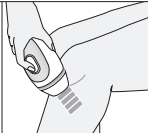
### Skin Test

Before your first treatment, we recommend you test your skin for reaction on a small area where you plan to treat. Perform two flashes next to each other as described below in chapter «Treatment methods», section «STAMPING».

Wait 15 minutes to examine the test area. If the skin appears normal, begin your first treatment. If a skin reaction develops, wait for 24 hours. If skin reaction is present after this period, discontinue using the device and consult your physician.

## Treatment methods

The area to be treated must be hair-free, clean and dry. Before pushing the treatment button, look away from the device to avoid the bright light.

	<p><b>GLIDING</b> for fast and efficient treatment of large areas (e.g. legs).</p>		<p><b>STAMPING</b> for treating small, hard-to reach areas (e.g. face for females, bikini line, underarms, knee, chin, ankles)</p>
<p>1. Place the treatment window (1) firmly on your skin. When the entire front area has skin contact, the power indicator (4) shows pink or white light to confirm that the device is ready to flash.</p>			
<p>2. <b>Press and hold</b> the treatment button (5). The device will begin to flash.</p>		<p>2. <b>Press and instantly release</b> the treatment button (5). The device will emit a flash.</p>	
<p>3. Move the device <b>slowly and continuously</b> from the bottom to the top of the treatment area. Make sure to always keep full skin contact so that the device flashes continuously.</p>		<p>3. Lift off the device and place it onto the next adjacent spot to be treated. Wait for the power indicator to light pink or white. Then press and release the treatment button.</p>	
<p>Once the gliding stroke is completed, lift the device off the treatment area. Repeat steps 1-3 to treat the next adjacent area in a parallel stroke.</p>		<p>Continue treatment spot by spot until you have completed the entire area you wish to treat.</p>	

Make sure you avoid missed areas by placing the treatment window directly next to the previously treated spot without leaving any space between treated areas. Avoid overlapping treatment or treating the same area twice as this may increase the chance of producing adverse reactions.

If you feel intense pain, stop treatment immediately and consult your physician before using the device again.

### How to treat sensitive areas

Before treating, remember to perform a skin test. Remove any visible hair and cosmetics, lotions or creams. You might want to confirm full hair removal using a mirror.

#### Face (females only)

When treating the upper lip area, avoid contact with the lips and the nostrils. The central area below the nostrils might be sensitive.

Do not use on permanent / semi-permanent lip liner on the lip margin.

TIP: You can press your lips together to ensure that the light does not activate on the upper lip itself.

Hair growth under the chin can be denser and coarser. This can sometimes lead to slight skin reddening following a treatment. This redness is temporary and should completely disappear after 24 hours.

#### Intimate area (refer to page 4 for suitable treatment areas)

Stand, sit or crouch down to reach all areas to be treated.

Apply the stamping method (spot-by-spot treatment).

**NOTE:**

- If redness or a burning sensation occurs and lasts for several minutes, apply a cool cloth/pad to soothe the area.
- If any pigment changes in the skin, e.g., lightening or darkening, are observed, stop using the device. Seek medical advice to determine the cause which may be due to an underlying medical condition.

## Treatment plan

Start-up phase	1 treatment per week (4-12 weeks)
Maintenance treatment	continue treating to your needs, e.g. every 1-2 months.

The typical hair growth cycle varies between individuals and could take up to 16 months depending on body area (underarm, lower leg, bikini area). Since only hairs in their growing phase are susceptible to treatment with light, it is important to do multiple treatments to achieve continuously smooth skin. For this reason, we recommend a start-up phase with one treatment per week over a period of 4 - 12 weeks.

For the facial area or those with darker skin tones, you may need to extend the start-up phase to achieve best results. In case you no longer see hair growing, you can switch to maintenance treatments.

## After treatment

Once you have completed your treatment, unplug the mains cord from the electrical outlet. After treatment, we recommend using sunscreen SPF 15 whenever treated areas may be exposed to the sun. We recommend avoiding exposure to artificial light sources, such as tanning booths or solariums.

Avoid any activity or skin products that may irritate your skin for 24 hours after treatment. This includes using hot tubs or saunas, wearing tight-fitting clothing, using bleaching creams, perfumed products or exfoliation products. The above may cause irritation, such as itching and redness, in the treated areas.

## Cleaning and storage

After each treatment, inspect the treatment window (1), the skin tone sensor (2), the glass filter (3) and the built-in reflector for damage.

Do not touch the glass filter immediately after use to allow it to cool down.

Wipe the treatment window with a dry, lint-free cloth if any marks or black spots are visible.

Do not rinse or immerse the device in water or any other liquid.

Store the device in a cool, dry place. Make sure the treatment window and skin sensor are protected from damage.

## Possible side effects

Some people will experience slight discomfort (e.g. heat or redness of the skin) when using the device – this is normal. You may feel heat and tingling when activating the flash, followed by redness in the skin. This type of reaction usually disappears within 24 hours. Adverse reactions are usually immediate or appear within 24 hours. In rare cases, they may take up to 72 hours to appear. If you experience abnormal skin reactions please contact your physician.

Side effect	How to assess and react
Warm feeling or tingling sensations during treatment, which typically disappear after a few seconds to a minute and decrease with continued use.	This is expected and normal for all IPL treatments. You can keep on using the device as instructed.
Discomfort or pain is intense during treatment or persists after a treatment.	Stop using the device and consult your physician before using it again.
Skin redness during or after treatment which disappears within minutes to several hours.	This is expected and normal for all IPL treatments. You can keep on using the device as instructed once skin redness has disappeared.
Skin redness which does not go away within 24 - 48 hours after treatment.	Stop using the device and consult your physician before using it again.

<b>In very rare cases:</b>	
Swelling and redness around the hair follicle that disappears within two or three days.	You can keep on using the device as instructed once swelling and redness around the hair follicle has disappeared.
Swelling and redness around the hair follicle that <b>does not</b> disappear within two to three days.	Stop using the device and consult your physician before using it again.
Temporary changes in skin color (lightening or darkening).	If your skin color changes, stop using the device and consult your physician.
Increase in hair growth in areas of treatment (more likely to be seen with subjects of Mediterranean descent).	If this increase persists over several treatments, stop using the device and consult your physician.

## Trouble shooting

The following troubleshooting guide may help you solving initial problems you may experience with this device.

Problem	Solution
Power indicator (4) does not light up when treatment window is placed on the skin.	Make sure the device is connected to an electrical outlet. Press the treatment button to activate it. Ensure that the entire treatment window area touches the skin. If needed, move the device to improve skin contact.
Flashing sequence slows down.	Device is in heat-protection mode. To get back to normal speed, allow device to cool down.
Red power indicator light	Invalid skin tone reading. Your skin is either too light or too dark. Device will not flash.
Red flashing indicator light	Device error. Device will not operate. Unplug and reconnect. If error persists, contact a Braun Service Centre for repair.
Power indicator light flashes alternately red and white.	End of product lifetime. Maximum number of flashes (300.000) has been reached.
The device is broken, cracked or appears defective.	<b>Do not</b> use if the handle, the base station or the mains cord is damaged. If you are in any doubt about the safety or suspect it is damaged in any way, you must not use it.

## Disposal

The product contains recyclable electronic waste. To protect the environment, do not dispose of it in the household waste, but take it to appropriate local collection points.

## Labels and symbols

The base station is marked with the following symbols:

	Type BF Applied Part		Serial number
	Class II equipment		WEEE waste disposal mark of the European Union

	Optical radiation		Conformity mark of the European Union
	Follow instructions for use		Manufacturer
	Not suitable for use in a bathtub or shower		Keep Dry

## Technical specifications

Repetition rate:	0.63 s, continuous operation
Max optical output:	4 J/cm <sup>2</sup>
Pulse length / profile:	0.5 - 1.0 ms FWHM
Wave length:	510 - 2000 nm
Power line input:	100 - 240V~, 50-60Hz, 1.7 - 0.91A
Operating temperature:	5°C - 40°C
Optimum operating temperature:	15°C - 30°C
Operating humidity:	up to 93% R.H. non-condensing
Operating pressures:	700 hPa to 1060 hPa
Treatment area (spot size):	3 cm <sup>2</sup>

### Transport Conditions:

Temperature range:	- 25°C to 70°C
Humidity:	up to 93% R.H. non-condensing
Pressures:	500 hPa to 1060 hPa

The device includes 300,000 flashes and is intended for use up to 5 years by a single user.

This device complies with the following EMC standards and latest amendments applicable to Household and Medical Devices. EN 55014-1:2017, EN 55014-2:2015, EN 61000-3-2:2014, EN 61000-3-3:2013, IEC 60601-1-2:2014 (Ed 4.0), EN 62311:2008, EN 62233:2008, FCC CRF 47 Part 15B, 15.107, 15.109, ICES-003 Issue 6.0.

This device also complies with the following safety standards and latest amendments applicable to Household and Medical devices: IEC 60335-1:2010/AMD2:2016, IEC 60335-2-113:2016, IEC 60601-1:2005/AMD1:2012, IEC 62304:2006, IEC 60601-1-11:2015, IEC 60601-1-6:2010/AMD1:2013, IEC 62366-1:2007/AMD1:2014, IEC 60601-2-57:2011, IEC 60601-2-83:2019, IEC 62471:2006.

This device is classified as Risk Group Exempt for all hazards.

For Spectral Irradiance, please refer to Annex 1 on page 97.

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