

## TECHNICAL NOTE

# Sequencing Metrics & Base Composition of Single Cell Multiome ATAC Libraries

## Introduction

The Chromium Next GEM Single Cell Multiome ATAC + Gene Expression workflow produces sequencing-ready Multiome ATAC and Gene Expression libraries from the same single nuclei, enabling simultaneous epigenomic and gene expression profiling. This Technical Note presents a comparison of sequencing metrics for a Single Cell Multiome ATAC library across multiple Illumina® platforms. Individual results may vary depending on the specific sequencing instrument and/or particular sample and loading characteristics.

As the Multiome Gene Expression library is identical to the Chromium Single Cell 3' Gene Expression (dual index) library, see the Gene Expression Dual Index Library data (Table 1) in the Technical Note—Sequencing Metrics & Base Composition of Single Cell 3' Gene Expression and Feature Barcode Dual Index Libraries (CG000374).

## Single Cell Multiome Libraries

Chromium Next GEM Single Cell Multiome ATAC + Gene Expression assay produces two library types from the same single nuclei:

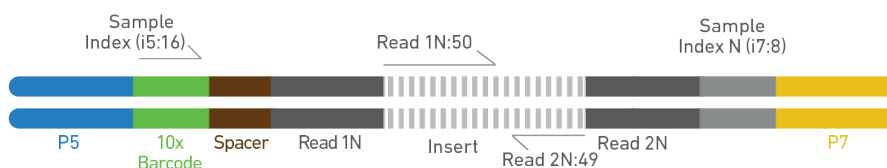
### 1. Chromium Single Cell Multiome ATAC Library

The library includes transposed DNA flanked by P5 and P7 sequences, necessary for binding to the Illumina flow cell. Read 1N and Read 2N sequences are used for priming and sequencing the DNA insert. The i7 index read is used for demultiplexing the libraries. The i5 index read is used to capture the 10x Barcode information.


### 2. Chromium Single Cell Multiome Gene Expression Library

The library includes cDNA insert flanked by P5 and P7 sequences, necessary for binding to the Illumina flow cell. TruSeq Read 1 is used to sequence the 16 bp 10x Barcodes and 12 bp UMI, and TruSeq Read 2 is used for priming and sequencing the cDNA insert. The two 10 bp sample indexes are sequenced in the i5 and i7 read respectively.

## Chromium Single Cell Multiome ATAC Library



## Chromium Single Cell Multiome Gene Expression Library

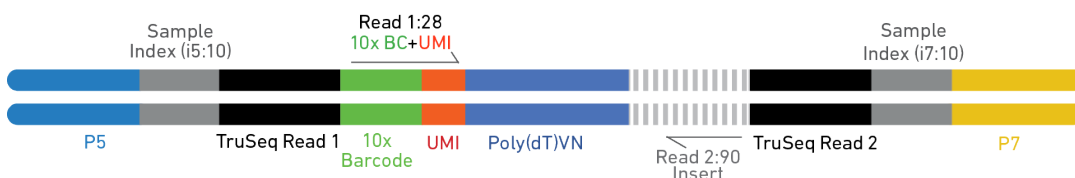

## Methods Overview

Multiome Single Cell ATAC and Gene Expression libraries were generated from one sample (PBMCs) as described in the Chromium Next GEM Single Cell Multiome ATAC + Gene Expression User Guide (CG000338). Approximately 500 nuclei were targeted using the Chromium Next GEM Chip J to generate the libraries used for this Technical Note.

As the Multiome Gene Expression library is identical to the Chromium Single Cell 3' Gene Expression (dual index) library, see the Gene Expression Dual Index Library data (Table 1) in the Technical Note—Sequencing Metrics & Base Composition of Single Cell 3' Gene Expression and Feature Barcode Dual Index Libraries (CG000374).

The Multiome ATAC library was quantified and sequenced with 1% PhiX on various Illumina® sequencers and their respective dual index workflows as listed below.

### 1. Forward Strand:

The i5 index Read occurs before Read 2 resynthesis, so the i5 index is sequenced on the forward strand. The sequencer will read the i5 index and capture the 10x Barcode in the Multiome ATAC library.

### 2. Reverse Complement:

The i5 index Read occurs after Read 2 resynthesis, which creates the reverse complement of the i5 index adapter sequence. Sequencers that run this workflow will read through the Spacer before reading through 10x Barcode in the Multiome ATAC library. This workflow requires a longer i5 index read length (24nt).

### 3. Custom Reverse Complement:


For sequencers that cannot support a longer i5 read of 24 cycles, a custom recipe is necessary. In this custom workflow, dark cycles skip sequencing the Spacer before reading the 10x Barcode in the Multiome ATAC library.

## Results Overview

Tables 1-3 show representative sequencing metrics and base composition data derived from Multiome ATAC library sequenced using the indicated sequencers/workflows. The Q30 quality scores, representative Data by Cycle plots, and other metrics for each sequencer/workflow is shown. Individual results may vary depending on the specific sequencing instrument and/or particular sample and loading characteristics.


### 1. Forward Strand Workflow

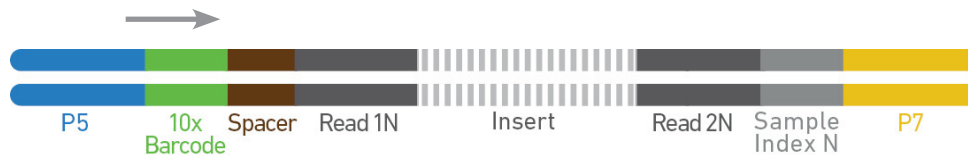Index primer binds P5 and reads towards Read 1N




### 2. Reverse Complement Workflow

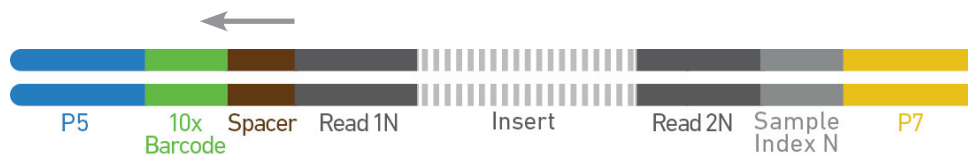Index primer binds Read 1N and reads towards P5



### 3. Custom Reverse Complement Workflow

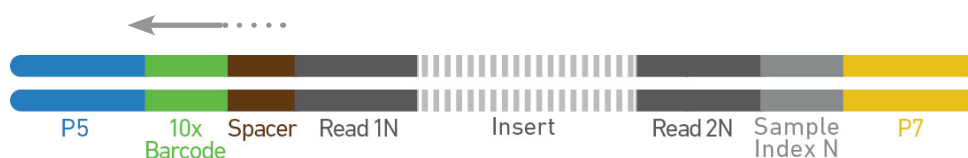Index primer binds Read 1N, performs 8 dark cycles (not reads), & sequences 10x BC



For more information about the workflows, visit the Illumina website and refer to the Indexed Sequencing Overview Guide (15057455 v07).

## Forward Strand Illumina® Workflow

Chromium Single Cell Multiome ATAC libraries were sequenced on Illumina sequencers with Forward Strand Workflow. 'Data by Cycle' plot from the Illumina software displaying the percentage of base calls and Q30 quality scores along with additional metrics are shown in Table 1.

Sequencing configuration & run parameters:



Paired-end, dual indexing




- Read 1: 50 cycles
- i7 Index: 8 cycles
- i5 Index: 24\* cycles
- Read 2: 49 cycles

\*Read length may be shortened to 16 bp

**Table 1: Representative Plots and Sequencing Data**

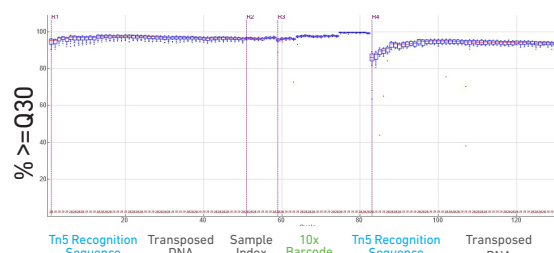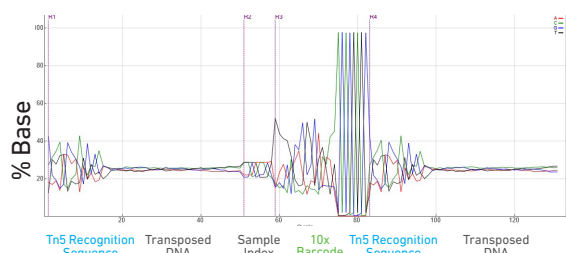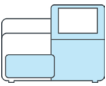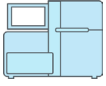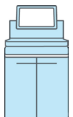Representative plots show the base distribution that reflects the Tn5 transposase recognition sequence, transposed DNA, sample index, 10x Barcode, and Spacer.



		% ≥Q30				Yield per Lane (Gb)		% ≥ mapQ30
		R1N	i7	i5	R 2N	R 1N	R 2N	
<b>MiSeq™</b>								
	Loading Conc. (pM): 10 Cluster Density: 900 K/mm <sup>2</sup> Phix (%): 1	97.72	97.61	97.00	95.39	0.99	0.97	86.69
<b>HiSeq™ 2500</b>								
	Loading Conc. (pM): 10 Cluster Density: 1,0416 K/mm <sup>2</sup> Phix (%): 1	97.35	96.28	95.84	94.30	8.67	8.49	85.70
<b>NovaSeq™ (S2 flow cell; v1 kits; v1.6 software)</b>								
	Loading Conc. (pM): 300 % PF: 79.70 Phix (%): 1	94.29	89.88	87.97	92.85	112.51	110.21	87.11

## Reverse Complement Illumina® Workflow

Chromium Single Cell Multiome ATAC libraries were sequenced on indicated Illumina® sequencer with Reverse Complement workflow. 'Data by Cycle' plot from the Illumina® software displaying the percentage of base calls and Q30 quality scores along with additional metrics are shown below.

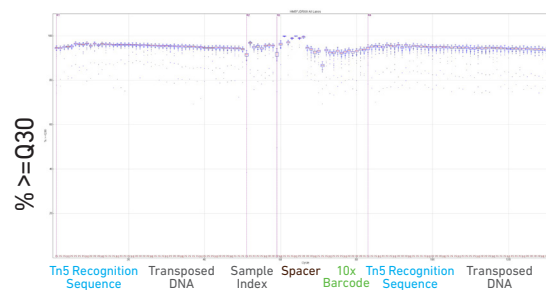
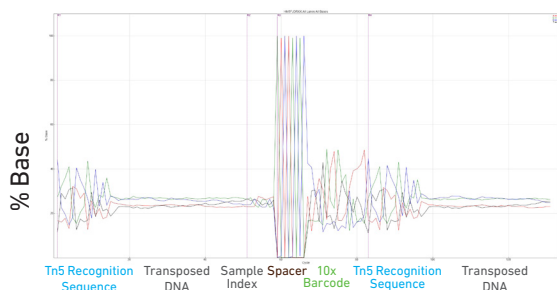
Sequencing configuration & run parameters:


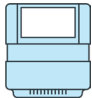
Paired-end, dual indexing

- Read 1: 50 cycles
- i7 Index: 8 cycles
- i5 Index: 24 cycles
- Read 2: 49 cycles

**Table 2: Representative Plots and Sequencing Data**

Representative plots show the base distribution that reflects the Tn5 transposase recognition sequence, transposed DNA, sample index, and 10x Barcode. A distinct 8 bp Spacer sequence is observed before reading the 10x Barcode.



		% ≥ Q30				Yield per Lane (Gb)		% ≥ mapQ30
		R 1N	i7	i5	R 2N	R 1N	R 2N	
<b>NovaSeq™ (SP flow cell; v1.5 kits; v1.7 software <i>released August 2020</i>)</b>								
	Loading Conc. (pM): 300 % PF: 72.27 Phix (%):1	95.42	94.84	94.43	94.68	22.63	22.16	85.37
<b>NextSeq™ 1000/2000*</b>								
	Loading Conc. (pM): 650 % PF: 83.03 Phix (%):1	94.91	90.43	90.66	93.11	34.49	33.78	87.94

\*NextSeq 1000/2000 data was derived from libraries that were generated separately from the libraries used in the other sequencers.

## Custom Reverse Complement Illumina® Workflow

Chromium Single Cell Multiome ATAC libraries were sequenced on the indicated Illumina® sequencer with Custom Reverse Complement workflow. 'Data by Cycle' plot from the Illumina software displaying the percentage of base calls and Q30 quality scores along with additional metrics are shown below.



Sequencing configuration & run parameters:


Paired-end, dual indexing

- Read 1: 50 cycles
- i7 Index: 8 cycles
- i5 Index: 8 **dark** cycles and 16 cycles
- Read 2: 49 cycles

**Table 3: Representative Plots and Sequencing Data**

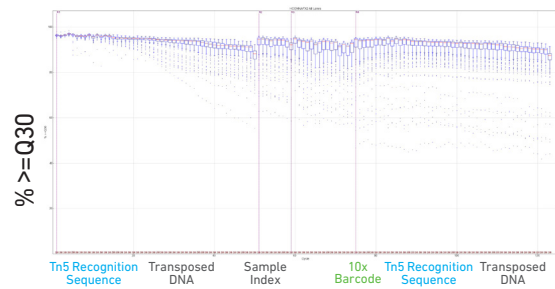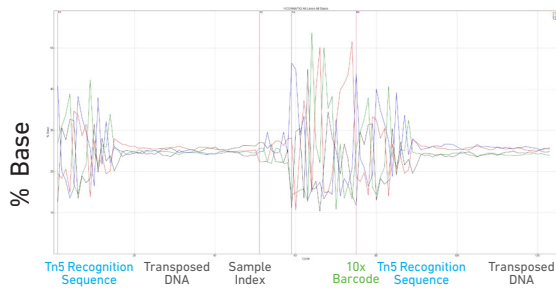Representative plots show the base distribution that reflects the Tn5 recognition sequence, transposed DNA, sample index and 10x barcode. A distinct 8 Spacer sequence is not observed prior to the 10x Barcode as dark cycles are not read by the sequencing instrument.



	% ≥Q30				Yield per Lane (Gb)		% ≥ mapQ30
	R 1N	i7	i5	R 2N	R 1N	R 2N	
<b>NextSeq™ 500/550</b>							
 Loading Conc. (pM): 1.5 Cluster Density: 205 K/mm <sup>2</sup> Phix (%):1	93.87	92.60	90.81	91.01	1.85	1.81	85.30

---

## Conclusions

In summary, 'Data by Cycle' plots displaying the percentage of base calls and Q30 quality scores% show highly consistent profiles for all sequencing platforms and workflows tested. The data serve as guidelines for assessing the quality of Single Cell Multiome ATAC library sequencing. Additional factors that may contribute to overall success of a sequencing run and impact downstream application performance metrics include:

- Sample preparation to obtain a high quality single nuclei suspension.
- Final libraries with fragment length of 200-2,000 bp and a significant number of inserts between 200-600 bp for optimal cluster formation on Illumina® flow cells.
- Reliable and accurate library quantification using the KAPA DNA Quantification Kit using the average insert size determined by Agilent Bioanalyzer QC.
- Sequencing platform loading concentration.

---

© 2021 10x Genomics, Inc. (10x Genomics). All rights reserved. Duplication and/or reproduction of all or any portion of this document without the express written consent of 10x Genomics, is strictly forbidden. Nothing contained herein shall constitute any warranty, express or implied, as to the performance of any products described herein. Any and all warranties applicable to any products are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. 10x Genomics provides no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein. The use of products described herein is subject to certain restrictions as set forth in the applicable terms and conditions of sale accompanying the purchase of such product. A non-exhaustive list of 10x Genomics' marks, many of which are registered in the United States and other countries can be viewed at: [www.10xgenomics.com/trademarks](http://www.10xgenomics.com/trademarks). 10x Genomics may refer to the products or services offered by other companies by their brand name or company name solely for clarity, and does not claim any rights in those third-party marks or names. 10x Genomics products may be covered by one or more of the patents as indicated at: [www.10xgenomics.com/patents](http://www.10xgenomics.com/patents). The use of products described herein is subject to 10x Genomics Terms and Conditions of Sale, available at [www.10xgenomics.com/legal-notices](http://www.10xgenomics.com/legal-notices), or such other terms that have been agreed to in writing between 10x Genomics and user. All products and services described herein are intended FOR RESEARCH USE ONLY and NOT FOR USE IN DIAGNOSTIC PROCEDURES.

The use of 10x Genomics products in practicing the methods set forth herein has not been validated by 10x Genomics, and such non-validated use is NOT COVERED BY 10X GENOMICS STANDARD WARRANTY, AND 10X GENOMICS HEREBY DISCLAIMS ANY AND ALL WARRANTIES FOR SUCH USE. Nothing in this document should be construed as altering, waiving or amending in any manner 10x Genomics terms and conditions of sale for the Chromium Controller or the Chromium Single Cell Controller, consumables or software, including without limitation such terms and conditions relating to certain use restrictions, limited license, warranty and limitation of liability, and nothing in this document shall be deemed to be Documentation, as that term is set forth in such terms and conditions of sale. Nothing in this document shall be construed as any representation by 10x Genomics that it currently or will at any time in the future offer or in any way support any application set forth herein.

Contact:  
[support@10xgenomics.com](mailto:support@10xgenomics.com)  
10x Genomics  
6230 Stoneridge Mall Road  
Pleasanton, CA 94588 USA

