

The new Lab Connect functionality enables you to independently include labs of your choice to your clinic. You can do this by browsing a list of iTero™-compatible labs registered in your scanner's country and appending the pertinent lab to the list accessible when completing the Rx.


Note: Available for users with administrator permissions only on MyiTero™ for doctors.

How to add a lab:

01

Select “Settings”

On the main page click the **Settings** button and select **Practice Settings**.

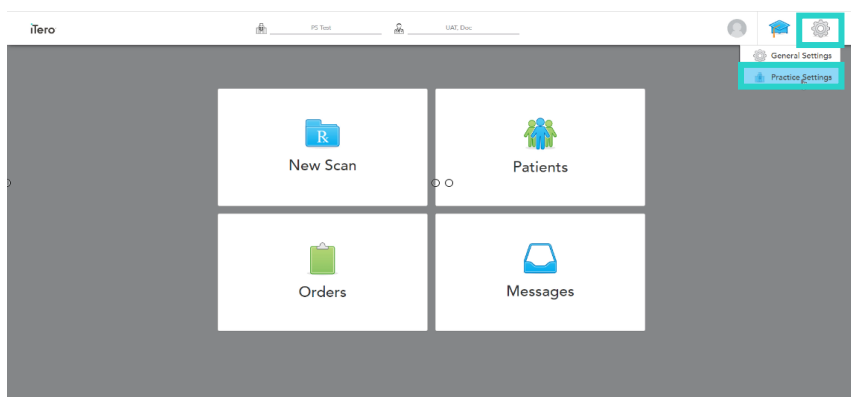


02

Select “My Labs”

Your list of Lab connections

- Click **My Labs**
- Click **Add Lab**


Note: If the **My Labs** tab is not displayed, contact your Admin to add the lab or to grant you admin permissions.

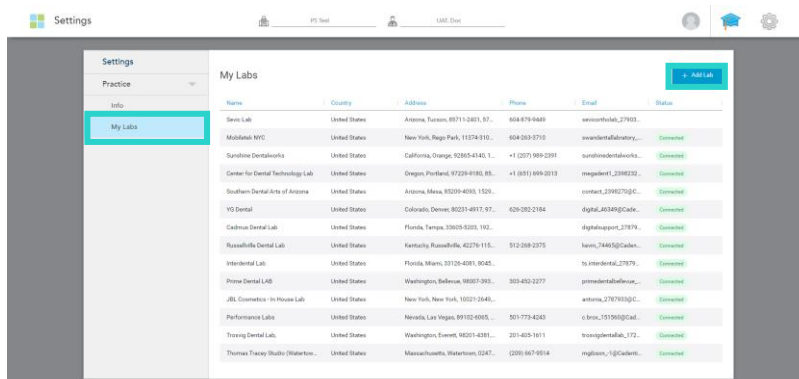

03

Find new lab


- Enter the name or part of a lab you want to work with and click **Search**.
- A list of labs that match your search criteria is displayed.

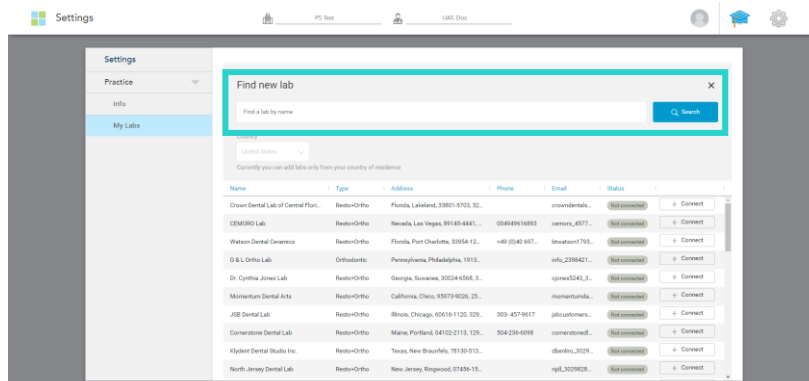

04

Connect to a new lab


- Click **Connect** to the lab you want to work with.

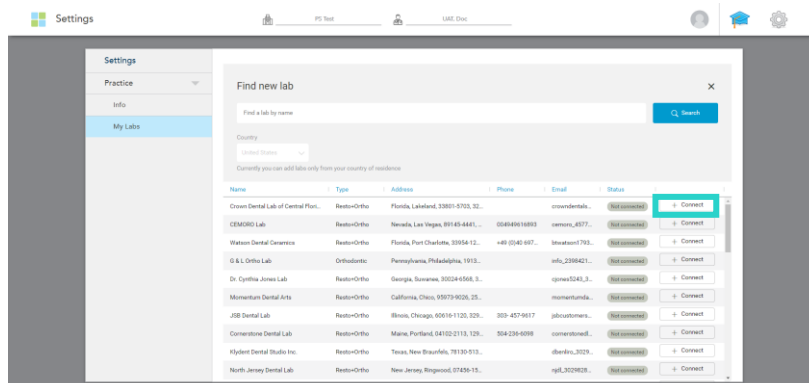

05

Confirmation


- A confirmation message is displayed, notifying you that the lab will be added to your account. The update will take a few minutes to appear in MyiTero™ portal.
- Click **Yes** to continue.

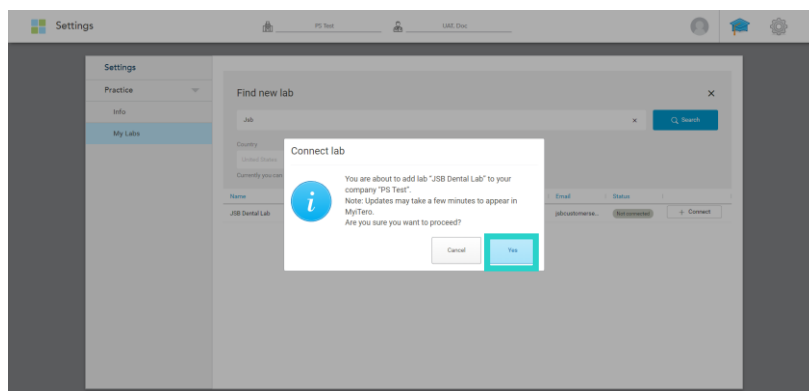

06

The Lab was added

- The recently added lab will appear at the top of the labs list in the **My Labs** tab. Its status will be **Pending** until the process is complete and it is linked to your account.

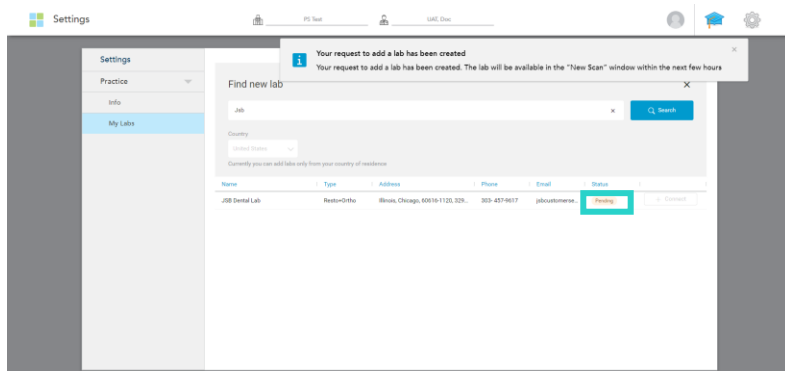

07

The Lab is connected

- Once the Lab is connected, it will appear under the **My Labs** tab with the status **Connected**.

