

Use of Safety Needle Devices Guidance

Medline Safety Syringe One Hand Activation

To enact safety feature, use one hand to gently press the safety shield against a hard surface.

HOW TO USE THE SAFETY SYRINGE WITH FLIP TECHNOLOGY



1. Remove the syringe from the packaging. Make sure the connection is secure by rotating the capped needle clockwise until snug.



2. Draw up and administer medication following the established protocol.



3. To enact the safety feature, use one hand to gently press the safety shield against a hard surface until you hear a click.




4. Place the syringe in the nearest sharps container for disposal.

Medline Safety Lancet: Blade retracts automatically



- Remove colored cap by twisting motion.
- Place lancet firmly on finger site. Push top lever bottom.
- Blade will advance and quickly retract into device.

Vacurette Safety Blood Collection Set

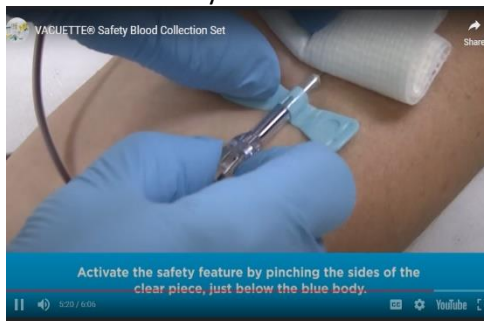

Click on Link for You Tube Video to watch how to activate the safety device by pinching the two buttons on either side device with one hand and pull to hear the click locking the sharp

[VACUETTE Safety Blood Collection Set - YouTube](#)

BD Vacutainer Push Button blood collection device for blood draws.
while the needle is still in the vein, activate the button with the tip of the index finger.
The needle will retract automatically from the vein. Apply pressure to the site.



BD Vacutainer Luer-Lok for Central line blood draws



Securely attach the BD Vacutainer Luer-Lok to the central line access port after disinfection.
After removing last tube, discard BD Vacutainer Luer-Lok device in approved sharps container.

One handed scoop if non-safety
Never a first choice

Place cap on table, hold used needle and syringe in one hand, then slide needle into cap. Lift-up and snap cap on.

