

Step-by-step guide: MyiTeroTM portal

Time saving features



- Fill in patient information or Rx remotely.
- Copy notes or Rx text from the computer.

Rx cannot be edited once the scanner is in scan mode and the order status is updated to scanning.

The screenshot shows the MyiTero portal interface. At the top, there's a header with 'Dentist: [Name]' and 'License: 12345'. Below this is a 'Patient:' section with fields for 'First Name:', 'Last Name:', 'Date of Birth:', 'Gender:' (with radio buttons for Male and Female), and 'Chart Number:'. There's a search icon and a 'Clear' button. Below the patient section is an 'Order:' section with 'Procedure *' (set to 'Fixed Restorative'), 'Type', 'Due Date:', and 'Send to:'. At the bottom is a 'Scan Options:' section with 'NIRI Capture:' (checked) and 'Pre-Treatment Scan:' (unchecked).

Educate patients

- Access patient's scans with the different views based on your need⁽¹⁾.
- Use the review tool⁽²⁾ with 5D imaging systems and/or iTeroTM NIRI (Near Infra-Red Imaging) technology and intraoral camera images.
- Use the snapshot tool⁽³⁾ to view and annotate⁽⁴⁾ prior scans.
- Annotated scans may be downloaded⁽⁵⁾, exported⁽⁶⁾, or viewed in the iTeroTM scan report.

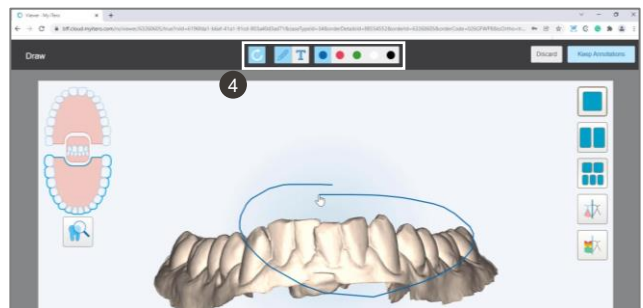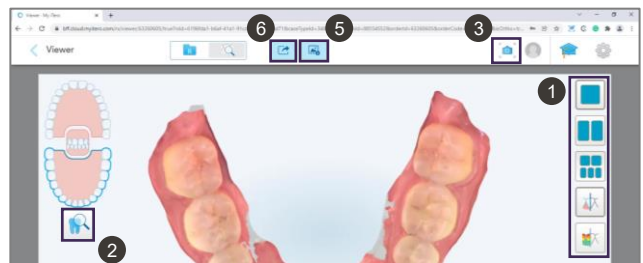


Exporting STLs, and downloading snapshots, and reviewing the scan reports are available only on the MyiTeroTM portal.

Step-by-step guide: MyiTero™ portal

Access programs

- Invisalign® Outcome Simulator, Invisalign® Progress Assessment, and OrthoCAD can be accessed from the MyiTero™ portal.
- Download the software on the computer.
- To open the Invisalign® Outcome Simulator, and the Invisalign® Progress Assessment applications, log in with your IDS (Invisalign® Doctor Site) to link the accounts.

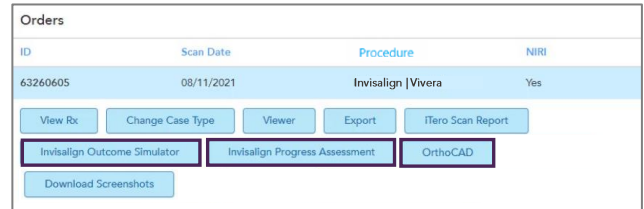

Account overview

- Change the password using Login.
- Set default Rx Settings.
- Adjust the Export Settings.
- Explore the scanner information under your account using System Settings.

