

**Thermo Clone™ VPS  
50 ml Cartridges**

50 ml	Regular 2pk	Fast 2pk
Super Light	—	4060
Light	4057	4065

2 x 50 ml cartridges  
12 x Each Mixing/IntraOral tips



50 ml	Regular 2pk	Fast 2pk
Medium	4058	4066
Heavy	4059	4067

2 x 50 ml cartridges  
6 x Mixing tips



50 ml	Fast 10pk
Light	4076

10 x 50 ml cartridges  
48 x Each Mixing/IntraOral tips



50 ml	Fast 10pk
Heavy	4077

10 x 50 ml cartridges  
30 x Mixing tips



**Thermo Clone™ Bite Registration**

50 ml	2pk
Bite Registration	4072

2 x 50 ml cartridges  
12 x Mixing tips



50 ml	2pk
Clear Bite Registration	4079

2 x 50 ml cartridges  
6 x Mixing tips



**Thermo Clone™ VPS  
380 ml Jumbo Cartridges**

380 ml	Regular 1pk	Fast 1pk
Medium	4068	4070
Heavy	4069	4071

1 x 380 ml cartridges  
10 x Dynamic Mixing tips



**Thermo Clone™ VPS Accessories**

**4073 - Thermo Clone VPS Putty Kit**  
1 x Each 200 ml base and 200 ml catalyst  
2 x Scoops



**5100 - Dispensing Gun 1pk**



**4074 - Tray Adhesive 10 ml Bottle 1pk**



Impression Mixing Tips	50pk
Yellow	2902
Teal	2904



IntraOral Impression Tips	50pk
	2906



Dynamic Mixing Tips	50pk
	4075



**TAKE A  
GREAT  
IMPRESSION  
THE FIRST  
TIME**




**THERMO  
CLONE™ VPS**

# DESIGNED TO CAPTURE EVERY DETAIL

Thermo Clone impression material is designed to capture even the finest features quickly and accurately. The results? Clear impressions with detailed margins—the first time. And because the temperature-sensitive material sets quickly, you'll see less distortion in your impressions and more smiles on your patients' faces.


## HYDROPHILIC

Hydrophilicity ensures precision in your impressions by displacing moisture, so the material can reach every part of every tooth. With a contact angle of 40 degrees, Thermo Clone material is among the industry leaders in hydrophilicity.

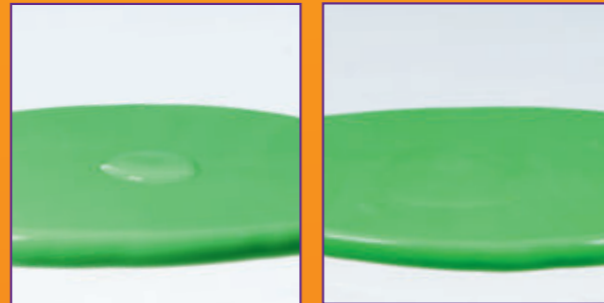

## THIXOTROPIC

Thixotropic materials become more fluid as they are agitated—like when they are applied to a crown preparation—and thicken when they are in place. This means that when Thermo Clone material is placed, the material flows into the sulcus and the spaces between teeth. This ensures a detailed impression and clear margins.

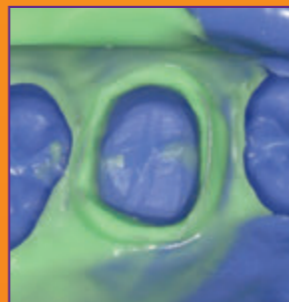

## EASY TO DELIVER

Thermo Clone impression material is available in 50 ml cartridges or 380 ml jumbo cartridges that fit most automatic mixing machines.

With an automatic mixing machine and the easy-to-use dynamic mixing tip, prepping to fill a tray is easy!



1. Remove cap      2. Attach Dynamic Mixing tip      3. Load into automatic mixing machine

SET TIMES	SET SPEED	WORKING TIME	SET TIME	TOTAL TIME
SUPER LIGHT BODY	FAST	1:00 min	1:15 mins	2:15 mins
LIGHT BODY	REGULAR	2:15 mins	2:15 mins	4:30 mins
	FAST	1:00 min	1:15 mins	2:15 mins
MEDIUM BODY	REGULAR	2:15 mins	2:30 mins	4:45 mins
	FAST	1:00 min	2:15 mins	3:15 mins
HEAVY BODY	REGULAR	2:15 mins	2:30 mins	4:45 mins
	FAST	1:00 min	2:15 mins	3:15 mins
PUTTY	FAST	2:00 min	2:00 mins	4:00 mins
BITE REGISTRATION	FAST	0:15 min	0:55 mins	1:10 mins

### Try it for yourself!


Apply Thermo Clone Light material to this square. When the material sets, lift it off the page to see the detailed impression.



# READY, SET, GO

Thermo Clone impression material is temperature sensitive. This means that as the temperature of the material increases, the setting time decreases. We call this a "Thermal-Accelerated Set".

At room temperature, Thermo Clone material stays malleable, with a working time of up to 1 min. Once the tray is placed in the patient's mouth, the material rapidly begins to set due to the increased temperature. This accelerated setting time means there's less chance of distortion. See below for a comparison of setting times for Thermo Clone heavy body fast set material.

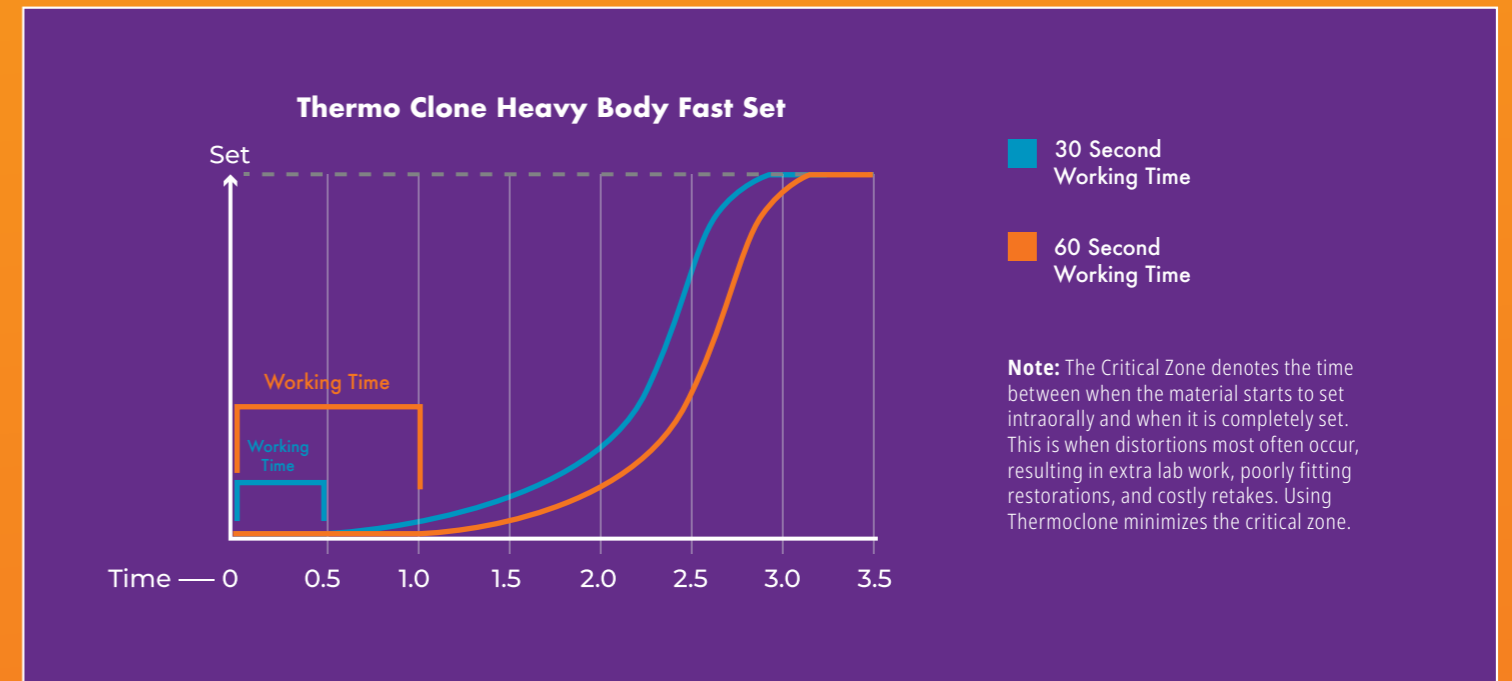

1. There is no minimum working time; Thermal-Accelerated Set ensures that the material begins to set as soon as it is placed in the patient's mouth.
2. In these examples, with 30 second and 60 second working times, ThermoClone was out of the critical zone and fully set at 2:45 mins and 3:15 mins respectively.

The results: a Thermal-Accelerated Set gives you a long working time, a short setting time, and minimal time in the critical zone.



**THERMO**  
**CLONE**™ VPS