

The *Kintrol* was designed for less active to moderately active users who mainly navigate indoor environments and use a walking aid outdoors, and place a high value on reliable support from their prosthetic foot.

Product features

- Hydraulic prosthetic foot
- Full length fibreglass spring provides the required stability and a soft rollover
- 12° range of motion (2° DF, 10° PF): increased range of motion for a more natural gait pattern and more stable standing
- The active dorsi-assist spring ensures ground clearance during the swing phase, which can help prevent the toes getting caught
- The Kintrol is lowered to the floor when sitting for a more natural look.
- Custom heel stiffness for a better balance thanks to four different heel wedge options
- Weatherproof

More than a foot. A foundation.

Technical data

K-Level / Mobility grade	1, 2
Max. body weight	300 lbs / 136 kg
Side	left (L), right (R)
Sizes	23-30 cm
Footshell colour	beige, light brown
Weight with footshell*	33.62 oz / 953 g
Heel height	3/8" / 10 mm
Build height*	4 7⁄16" / 112 mm

* Reference size is 26 cm.

ottobock.

Easy rollover at slow walking speeds

Increased ground contac improves balance

Adjustable heel

ottobock.

VS4 Order form.


Customer

Customer no.	
Company	
Company	
Street	
Street	
Postal code	City
Postal code	City


O&P professional

Selection of the spring stiffness relative to


Body weight (kg)	Stiffness
45 – 59	1
60 - 81	2
82 – 109	3
110 – 136	4

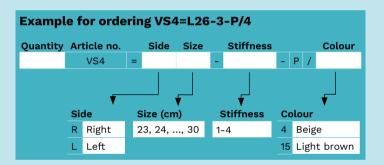

Alternative ordering option - Article number foot module (without footshell)

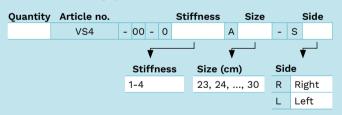

Alternative ordering option - Article number footshell

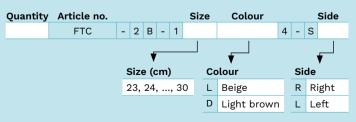


Compatible mechatronic knee joints

Kenevo

Shipping address (if different from customer address)

Customer no.	
Company	
Street	
Postal code	City
Customer ref.	

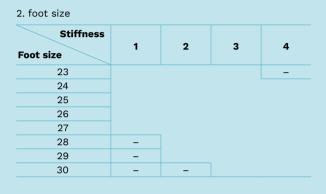

Scope of delivery

The scope of delivery includes the VS4* foot module, the FTC-2B-1* footshell, the 2F13=* heel wedge set (comprising four different degrees of hardness) and a black Spectra sock. Additionally, a tool for changing the footshell can be ordered (ACC-00-10300-00).

