

Genu Arexa PCL. Always the right support.

The **Genu Arexa PCL** knee brace is a device that is used for acute and follow-up treatment after a posterior cruciate ligament tear. The brace can support the healing process after an operation (ligament reconstruction) or as part of conservative therapy without an operation. The **Genu Arexa PCL** is also used after distensions or partial ruptures of the posterior cruciate ligament.

Indications

Rupture of the posterior cruciate ligament (acute treatment, conservative therapy, postoperative treatment)

Posterior cruciate ligament tear

Generally, a posterior cruciate ligament tear is caused by the application of a powerful external force, for example, due to accidents, contact sports or severe falls. Tearing of the ligament leads to recurring posterior displacement of the lower leg, which destabilises the leg.

Two variants for customised patient treatment

The **Genu Arexa** is a device for the conservative treatment of cruciate ligament injuries such as distensions, partial ruptures or ruptures.

The **Genu Arexa PCL** is especially intended for treating the posterior cruciate ligament.



50K13
Genu Arexa




50K14
Genu Arexa PCL



50K14

Genu Arexa PCL.

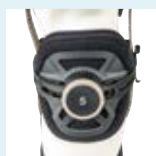

Provides stability after a posterior cruciate ligament tear

The brace was specially designed to treat injuries or tears of the posterior cruciate ligament.



Dynamic Pad

The large calf pad dynamically adapts to the leg's position thanks to an integrated spring.



Click-2-Go system

The Click-2-Go system makes it possible to adjust movement limits quickly and without tools. If the knee is swollen, the brace can be adjusted to the width of the individual's knee using spacers.




Easy to put on

The ***Genu Arexa PCL*** is applied over the front side of the leg. This makes it easier to put on.



Auto-adaptive plastic wings

The plastic wings on the top and bottom ends of the brace are adjusted to the contour of the individual's leg with the straps, making the brace more comfortable to wear.

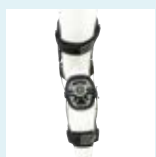

Ready to fit

The straps and pads of the ***Genu Arexa PCL*** are pre-adjusted to facilitate the initial adaptation of the brace.

Order information:

Art. no. 50K14 | Order no: article number=side-size | Example for ordering: 50K14=R-XXL

Side	Size	Circumference A (cm)	Circumference B (cm)
L/R	S	32–37	39–46
L/R	M	37–41	46–53
L/R	L	41–46	53–61
L/R	XL	46–50	61–69
L/R	XXL	50–56	69–79

