How to Deliver a Bolus


Open Printable Version

Note: Images are for educational purposes only. Your screen may look different.




Not all steps are shown


Follow the additional in-app prompts to deliver a bolus.

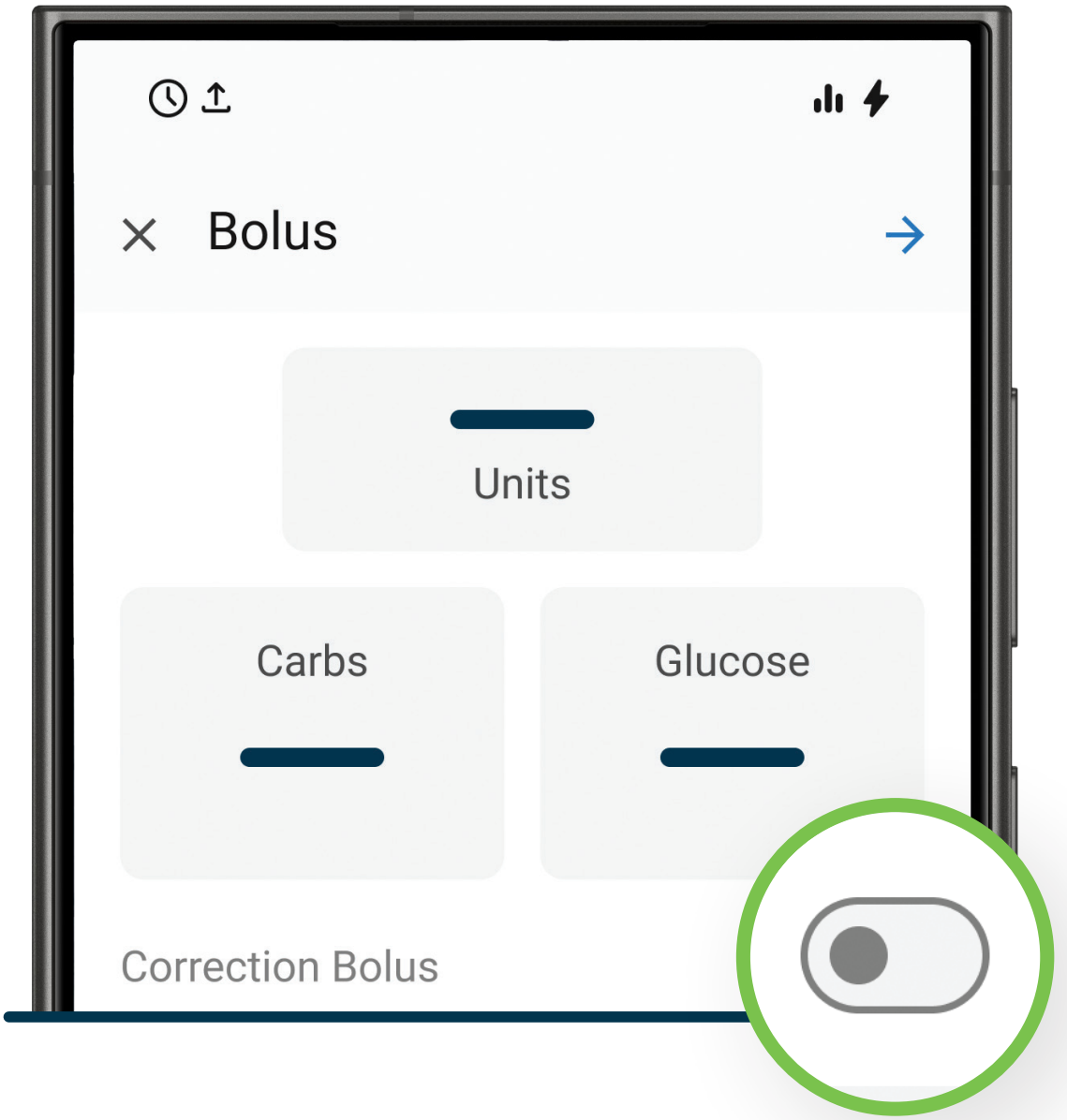


1 Use the toggle to turn the **Correction Bolus** off or on.

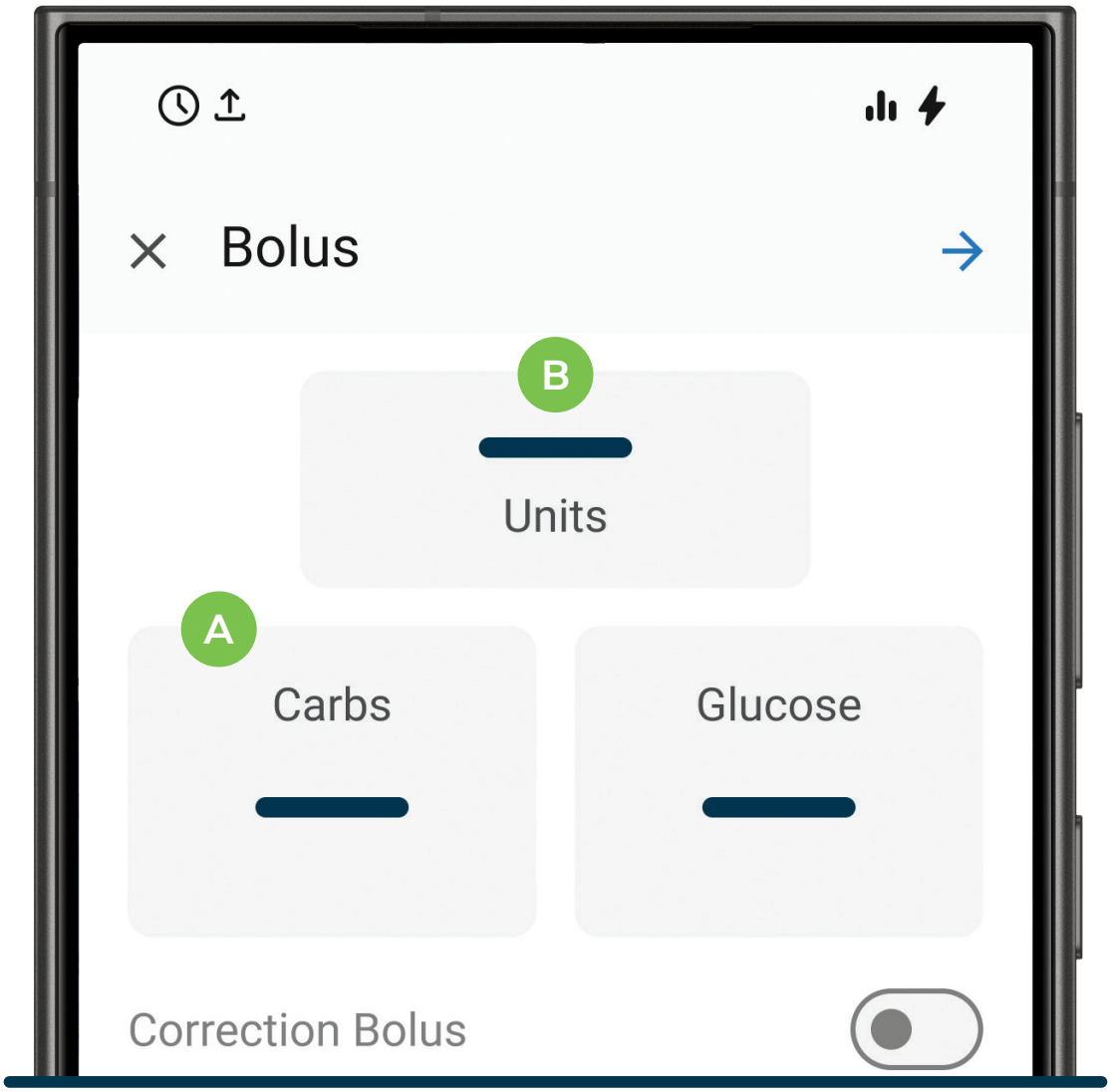


❷ Enter **Carbs** Ⓐ to calculate your dose, or manually enter the **Units** Ⓑ to deliver.

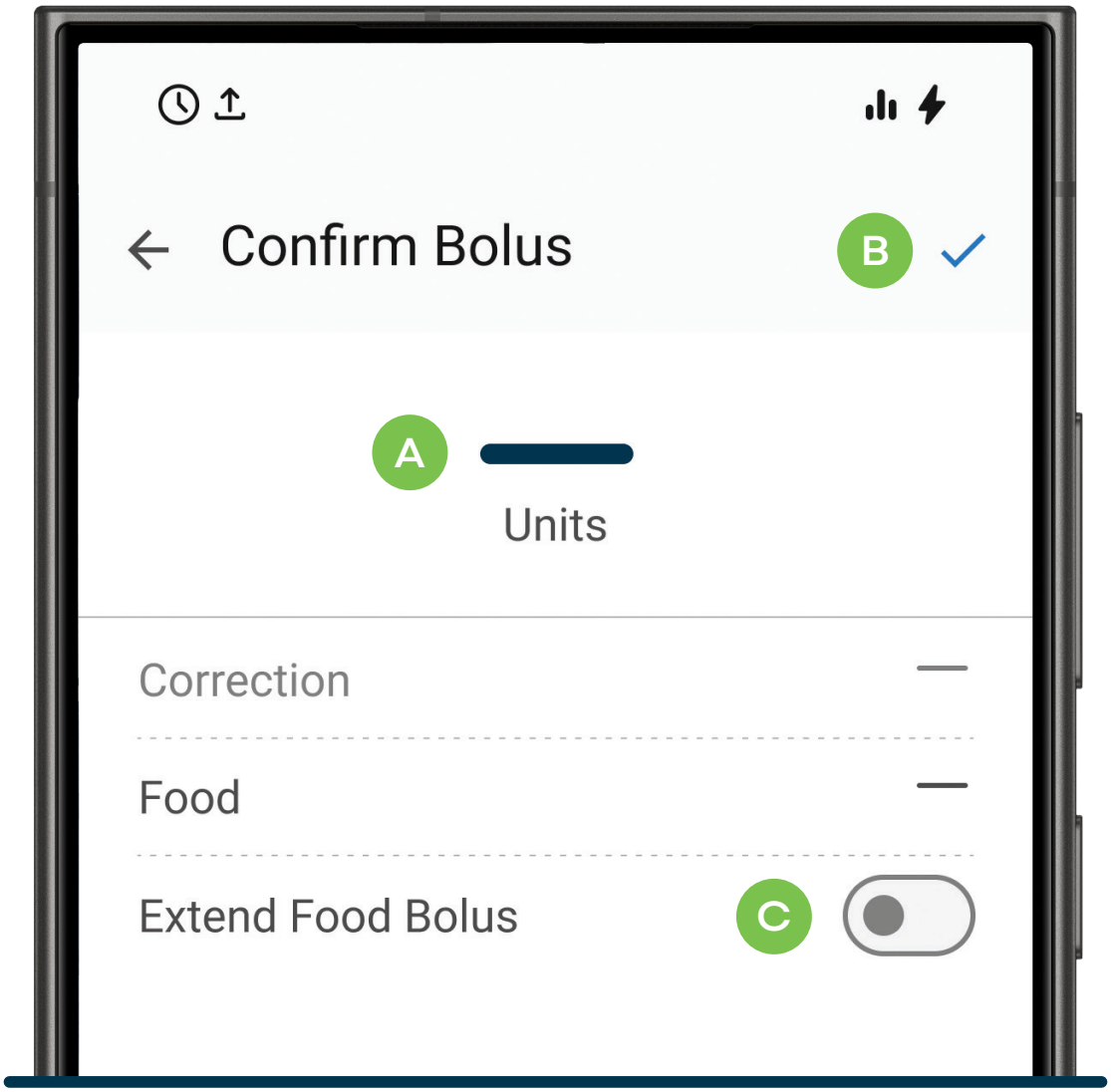

3 Confirm the unit dose A and follow the prompts B to complete the delivery. Or use the toggle C to extend the bolus.

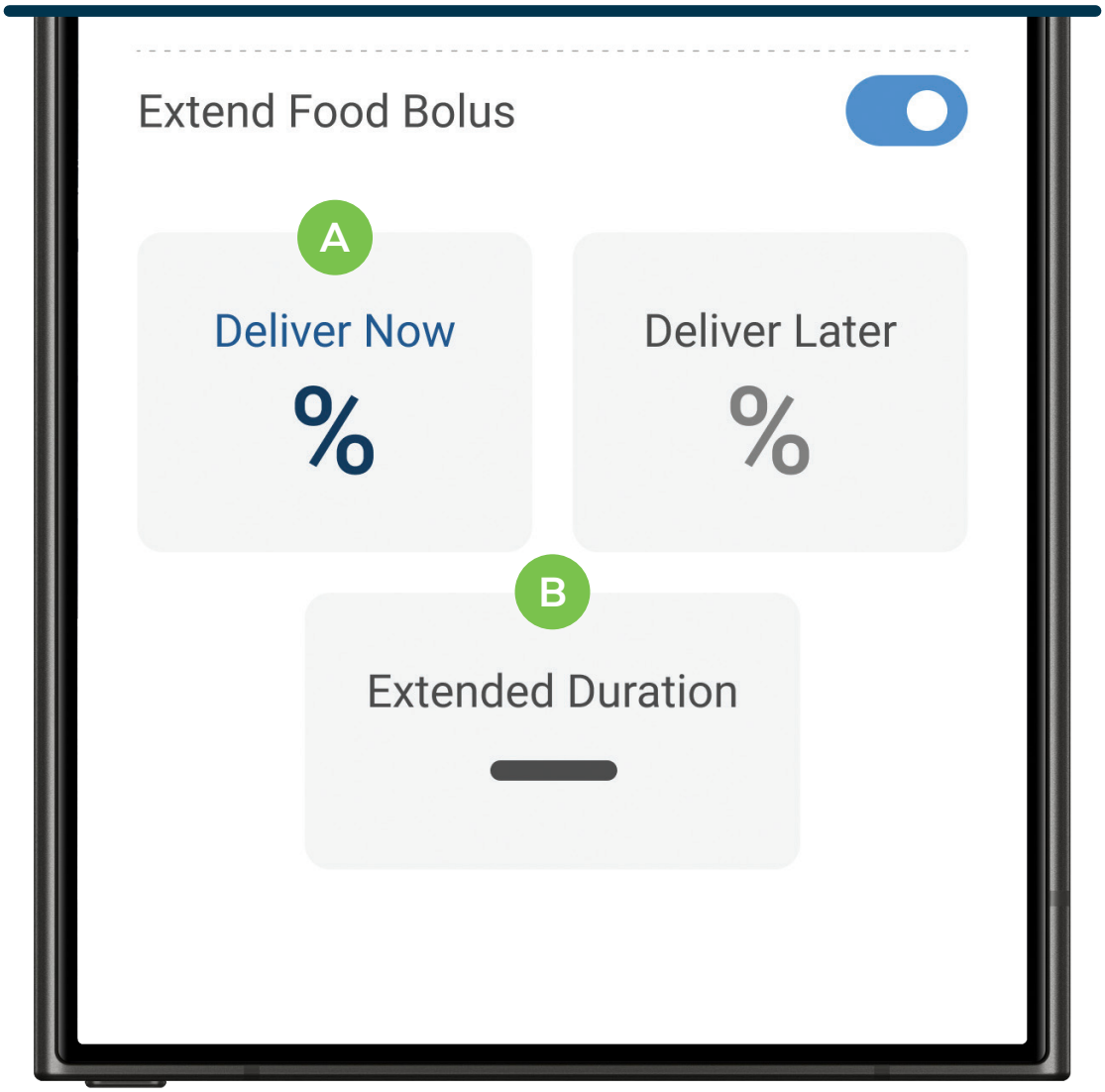

4 Set how much to **Deliver Now** A and the **Extended Duration** B for the rest.

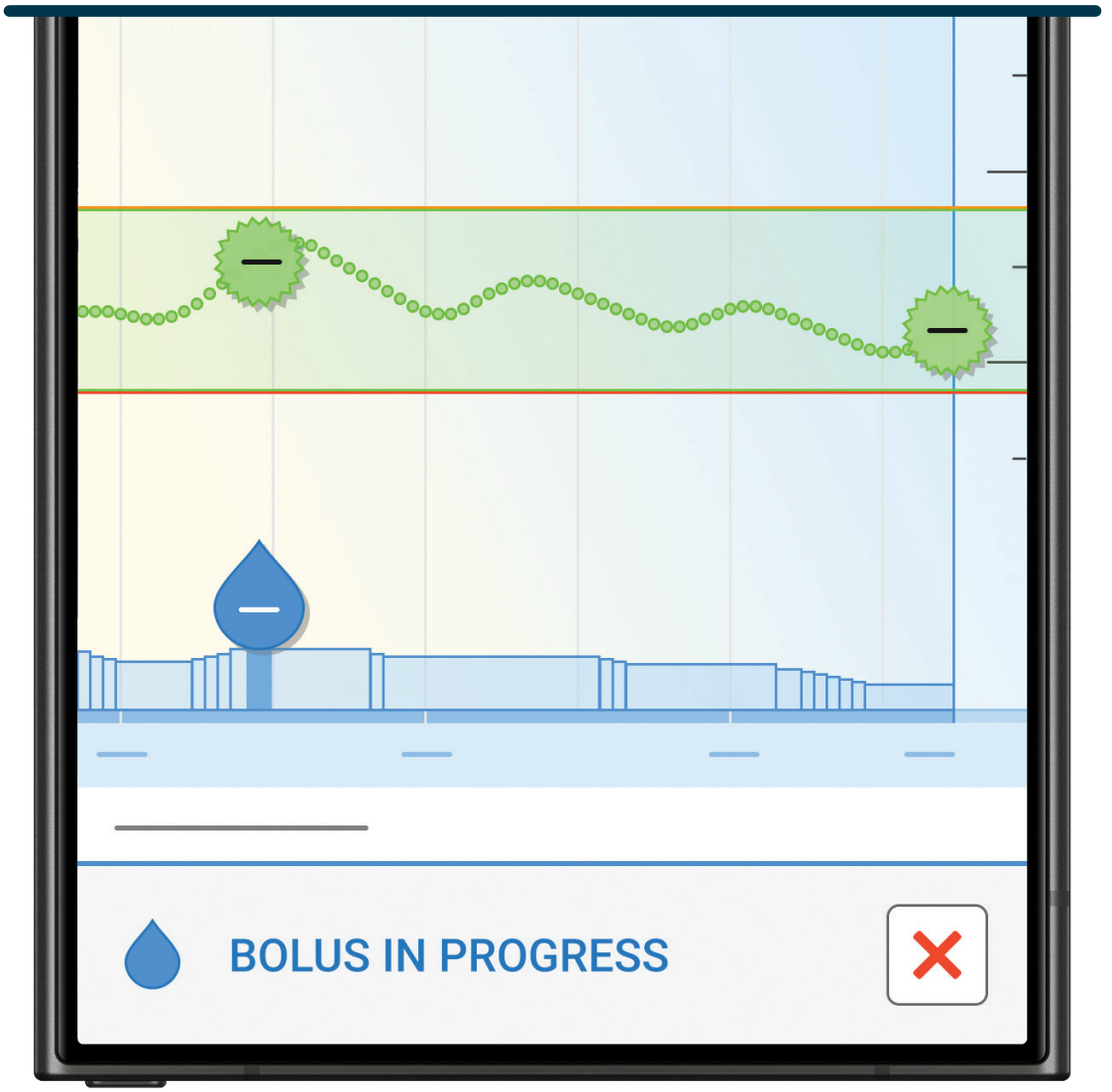

Tap the red “X” on the Dashboard to cancel at any time. The graph and IOB will reflect any insulin delivered.

These instructions are provided as a reference tool for individuals who are already familiar with the use of a Tandem Diabetes Care insulin pump. For more detailed information on the operation of your system, please refer to its user guide.

Important Safety Information

RX ONLY. Tandem insulin pumps are intended for the subcutaneous delivery of insulin, at set and variable rates, for the management of diabetes mellitus in persons requiring insulin.

Tandem insulin pumps with Control-IQ+ technology are intended for the management of type 1 diabetes in individuals 2 years of age and greater and of type 2 diabetes in persons 18 years of age and greater. Control-IQ+ technology is intended for use in pregnancy complicated by type 1 diabetes mellitus, provided the linked CGM system is suitable for use in pregnancy. Visit tandemdiabetes.com/safetyinfo for additional important safety information.

© 2026 Tandem Diabetes Care, Inc. All rights reserved. Tandem Diabetes Care, Tandem logos, and Control-IQ+ are either registered trademarks or trademarks of Tandem Diabetes Care, Inc. in the United States and/or other countries. MAT-1204_2

TANDEM
Diabetes Care