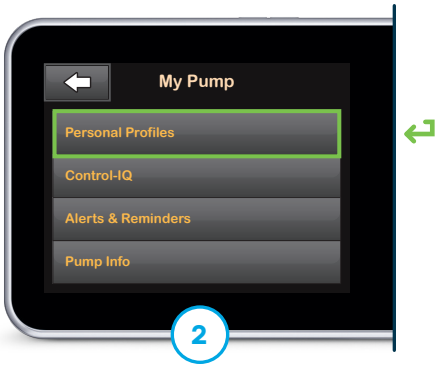
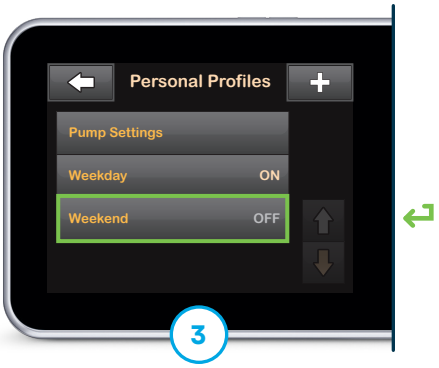


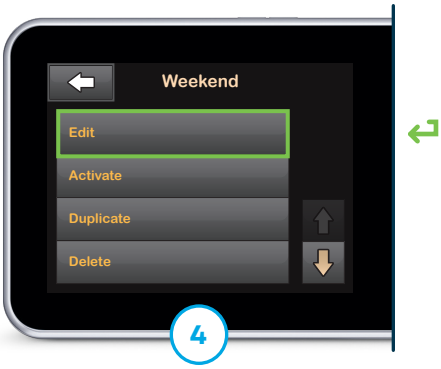

Tap **Options** from the Home screen and then **My Pump**




Tap **Personal Profiles**




Tap the profile name that needs changes




Tap **Edit** to make changes



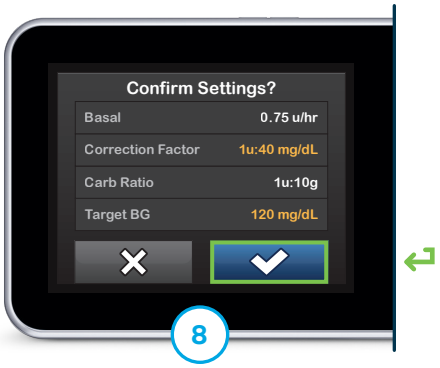

Tap **Timed Settings**



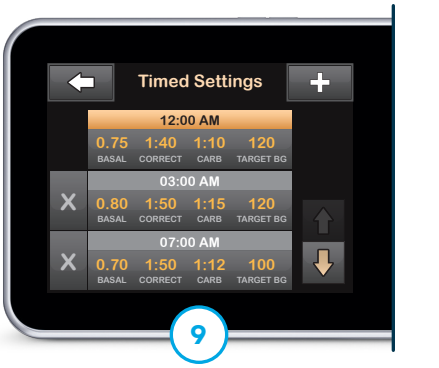

Tap the desired segment and then the setting that that needs to be changed




Use the keypad to enter the new value. Tap ✓ to save




Tap ✓ to save the segment and then tap ✓ to confirm

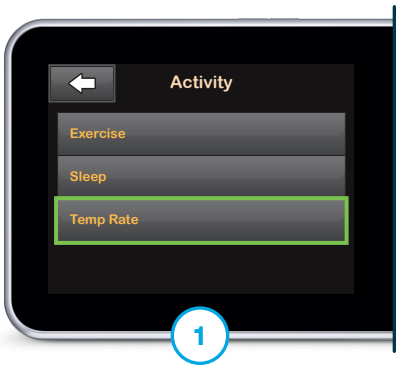



Tap other segments to make changes or return to the Home screen

# Program a Temp Rate



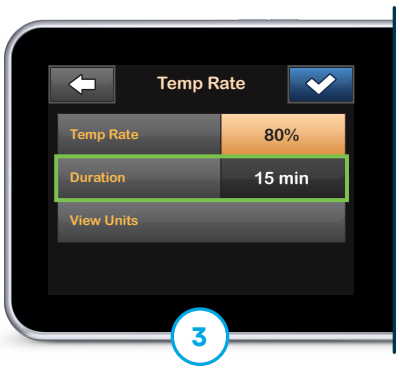
A Temp Rate is used to increase or decrease the current Personal Profile's Basal Rate for a period of time.




Tap **Options** from the Home screen, **Activity,** and then **Temp Rate**




Tap **Temp Rate,** enter the desired percentage, and then tap ✓ to continue




Tap **Duration,** enter the desired length of time, and then tap ✓ to continue. Tap ✓ again to confirm.




While active, a red or orange icon will display on the Home screen




To stop at any time, tap ✕ on the **Activity** screen and then ✓ to confirm

**Visit tandemdiabetes.com/support for additional resources.**

**These instructions are provided as a reference tool for individuals who are already familiar with the use of a Tandem Diabetes Care insulin pump. Not all screens are shown. For more detailed information on the operation of the t:slim X2 insulin pump, please refer to its user guide.**

**Important Safety Information:** RX ONLY. Indicated for patients with type 1 diabetes, 2 years and older and for patients with type 2 diabetes, 18 years and older. Safety info: tandemdiabetes.com/safetyinfo.

**WARNING:** Control-IQ+ technology should not be used in anyone under the age of 2 years old with type 1 diabetes or under the age of 18 years old with type 2 diabetes. It should also not be used in patients who require less than a total daily insulin dose of 5 units of insulin per day or who weigh less than 20 pounds, as those are the required minimum values needed for Control-IQ+ to operate safely.

© 2025 Tandem Diabetes Care, Inc. All rights reserved. Tandem Diabetes Care, Tandem logos, t:slim X2, and Control-IQ+ are either registered trademarks or trademarks of Tandem Diabetes Care, Inc. in the United States and/or other countries. MAT-0199_1