How to View TDI

Open Printable Version

Note: Images are for educational purposes only. Your screen may look different.



Not all steps are shown


Use this document as a guide and complete any additional steps as needed.



Viewing total daily insulin (TDI) on the Tandem Source platform (source.tandemdiabetes.com) is available for both the Tandem Mobi system and t:slim X2 insulin pump.


Direct viewing on the t:slim X2 insulin pump is also available.

TANDEM SOURCE

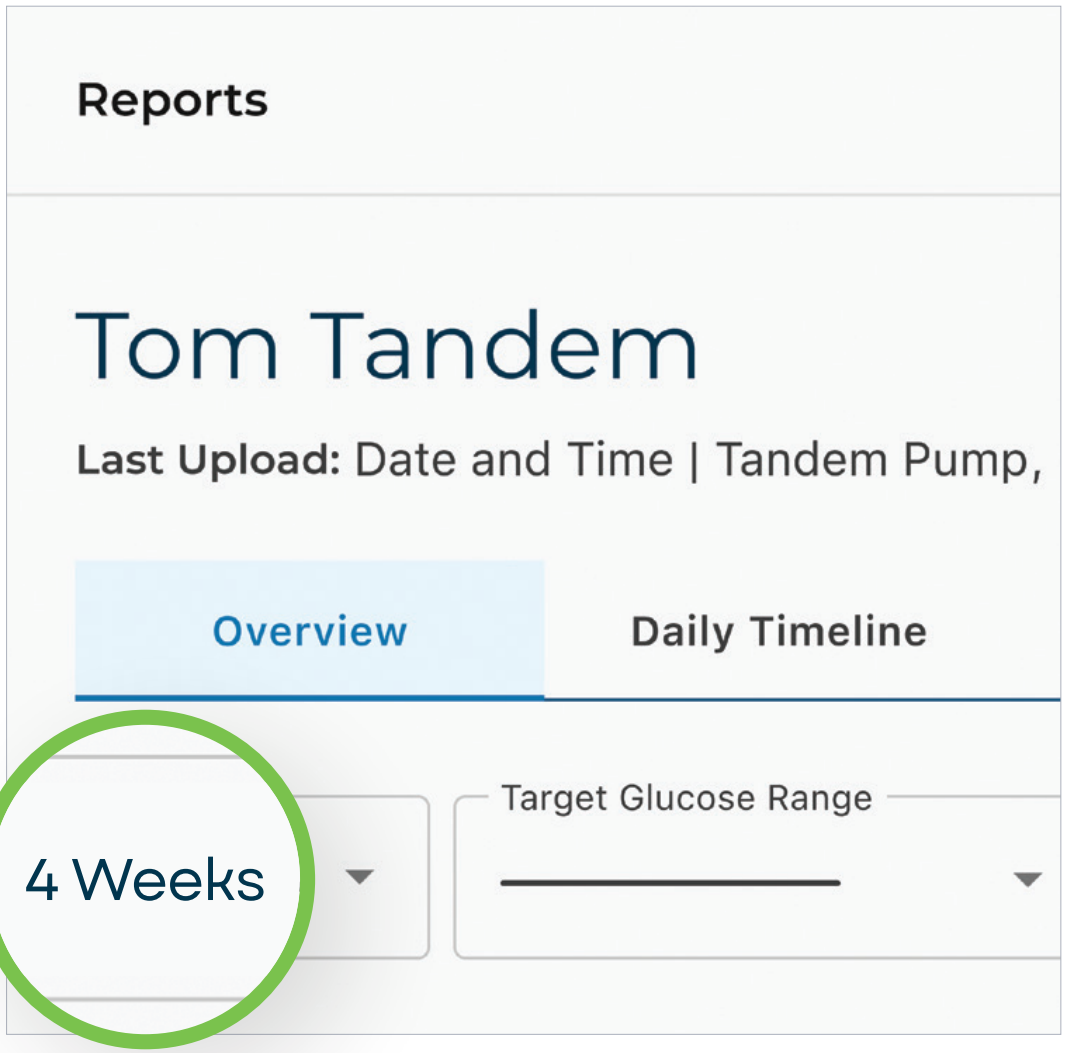


❶ Click on the **Reports** icon. Select 4 Weeks on the Overview tab to view a monthly summary.

TANDEM SOURCE

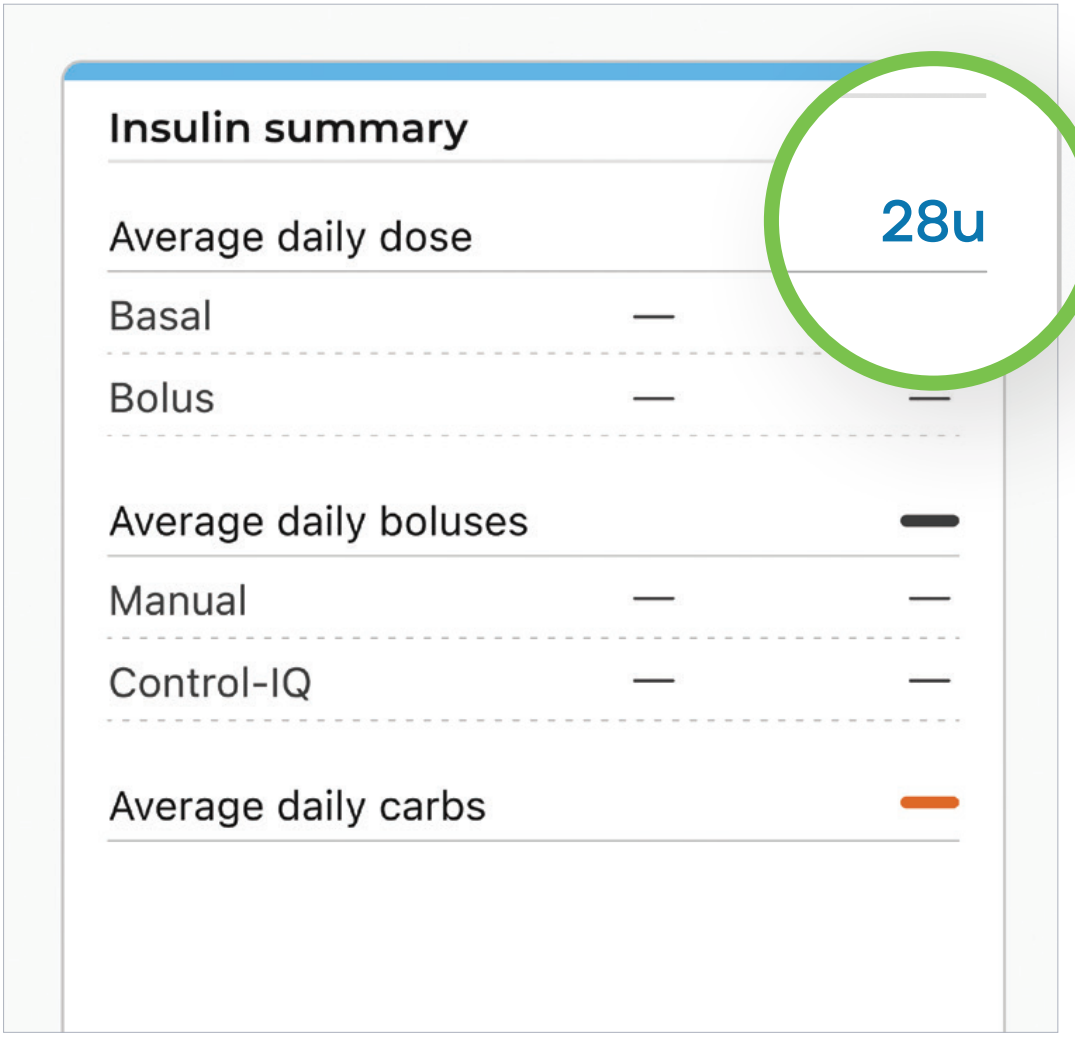

❷ Scroll to **Average daily dose** (average TDI for the selected range) on the Insulin Summary card.

TANDEM SOURCE

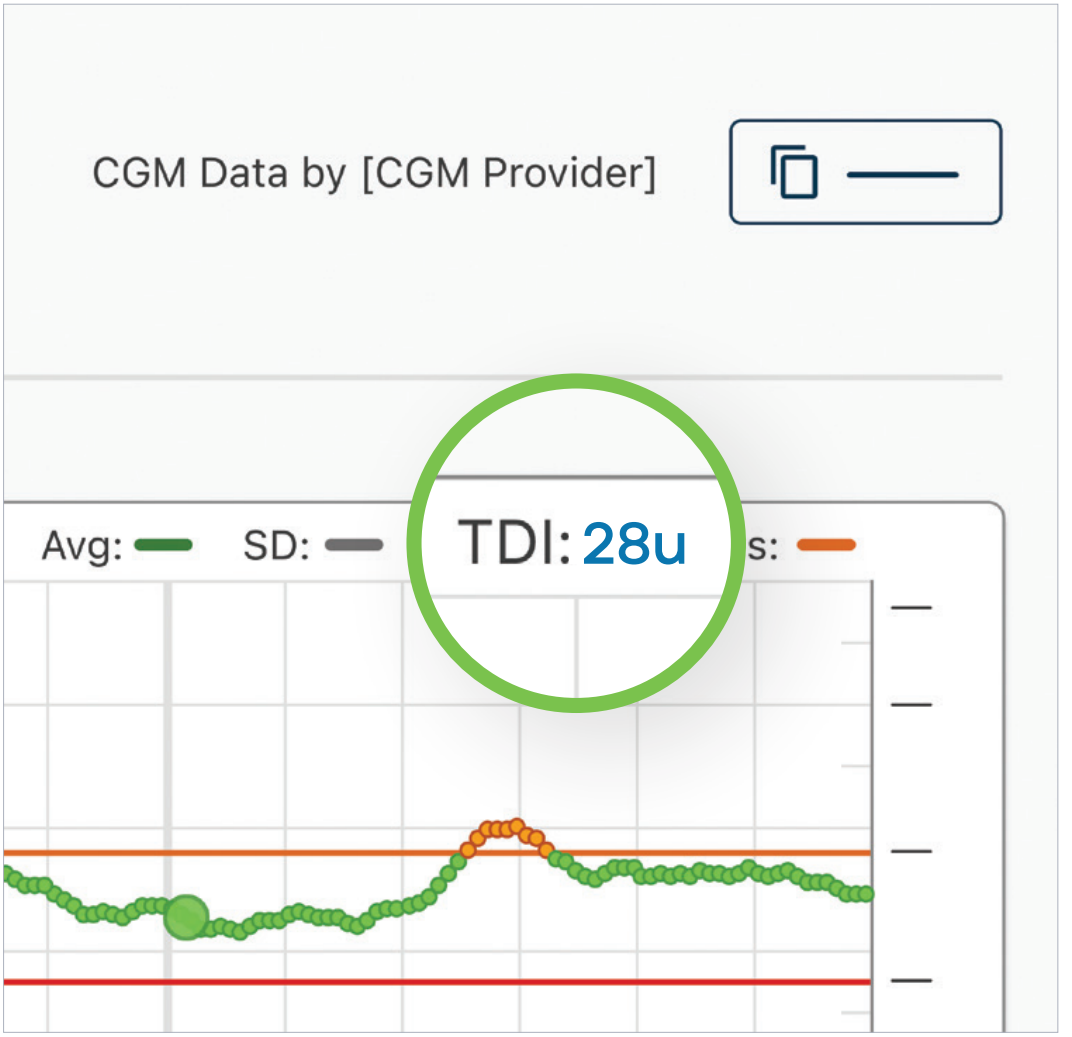

3 Select the **Daily Timeline** report tab to view TDI for a specific day.

T:SLIM X2



Tap Options, History, Pump History, then **Delivery Summary** (date range) or **Total Daily Dose** (specific day).

These instructions are provided as a reference tool for individuals who are already familiar with the use of a Tandem Diabetes Care insulin pump. For more detailed information on the operation of your system, please refer to its user guide.

Important Safety Information

RX ONLY. Tandem insulin pumps are intended for the subcutaneous delivery of insulin, at set and variable rates, for the management of diabetes mellitus in persons requiring insulin.

Tandem insulin pumps with Control-IQ+ technology are intended for the management of type 1 diabetes in individuals 2 years of age and greater and of type 2 diabetes in persons 18 years of age and greater. Control-IQ+ technology is intended for use in pregnancy complicated by type 1 diabetes mellitus, provided the linked CGM system is suitable for use in pregnancy. Visit tandemdiabetes.com/safetyinfo for additional important safety information.

The Tandem Source platform is intended for use by individuals with diabetes mellitus who use Tandem Diabetes Care insulin pumps, their caregivers, and their healthcare providers in home and clinical settings. The Tandem Source platform supports diabetes management through the display and analysis of information uploaded from Tandem insulin pumps.

© 2026 Tandem Diabetes Care, Inc. All rights reserved. Tandem Diabetes Care, Tandem logos, Tandem Source, Tandem Mobi, t:slim X2, and Control-IQ+ are either registered trademarks or trademarks of Tandem Diabetes Care, Inc. in the United States and/or other countries. MAT-2615_1

TANDEM
Diabetes Care