



**4018-1 - Jiffy Universal Extraoral Adjusting & Polishing Kit**  
 2 x Jiffy HP Grinder Wheels (1 coarse, 1 medium)  
 2 x Jiffy HP Grinder Tapers (1 coarse, 1 medium)  
 2 x Jiffy HP Universal Wheels (1 medium, 1 fine)



**4019-1 - Jiffy Universal Intraoral Adjusting & Polishing Kit**  
 2 x Jiffy RA Universal Points (1 medium, 1 fine)  
 2 x Jiffy RA Universal Cups (1 medium, 1 fine)  
 2 x Jiffy RA Natural Universal Wheels (1 medium, 1 fine)



**4236-1 - Jiffy HP Medium Polishing Wheel 1pk**  
**4237-1 - Jiffy HP Fine Polishing Wheel 1pk**



**4108-1 - Jiffy RA Medium Polishing Point 5pk**  
**4109-1 - Jiffy RA Fine Polishing Point 5pk**



**4238-1 - Jiffy HP Coarse Polishing Grinder Wheel 1pk**  
**4241-1 - Jiffy HP Medium Polishing Grinder Wheel 1pk**



**4234-1 - Jiffy RA Medium Polishing Cup 5pk**  
**4235-1 - Jiffy RA Fine Polishing Cup 5pk**



**4239-1 - Jiffy HP Coarse Polishing Grinder Taper 1pk**  
**4242-1 - Jiffy HP Medium Polishing Grinder Taper 1pk**



**1029 - Jiffy Goat Hair Brush 10pk**



**5540 - Diamond Polish Mint Syringe 0.5 µm 2pk**  
**5541 - Diamond Polish Mint Syringe 1 µm 2pk**  
 1.2 ml syringes

For a final esthetic finish, try the unique Jiffy Goat Hair Brush and Ultradent Diamond Polish Mint.


For more information about Ultradent's Jiffy™ Universal Ceramic finishing systems, please contact your sales representative!



SCAN TO VIEW PRICING

## UNIVERSAL. ESSENTIAL. EFFICIENT.

Ultradent's Jiffy Universal Ceramic Polishing System is a finishing solution in which all adjusters and polishers work interchangeably on porcelain, lithium disilicate, and zirconia. Each convenient kit includes everything you need, eliminating the need to purchase unnecessary or multiple components and allowing you to create a truly smooth, beautiful, high-gloss finish of the highest quality on all your ceramic restorations.

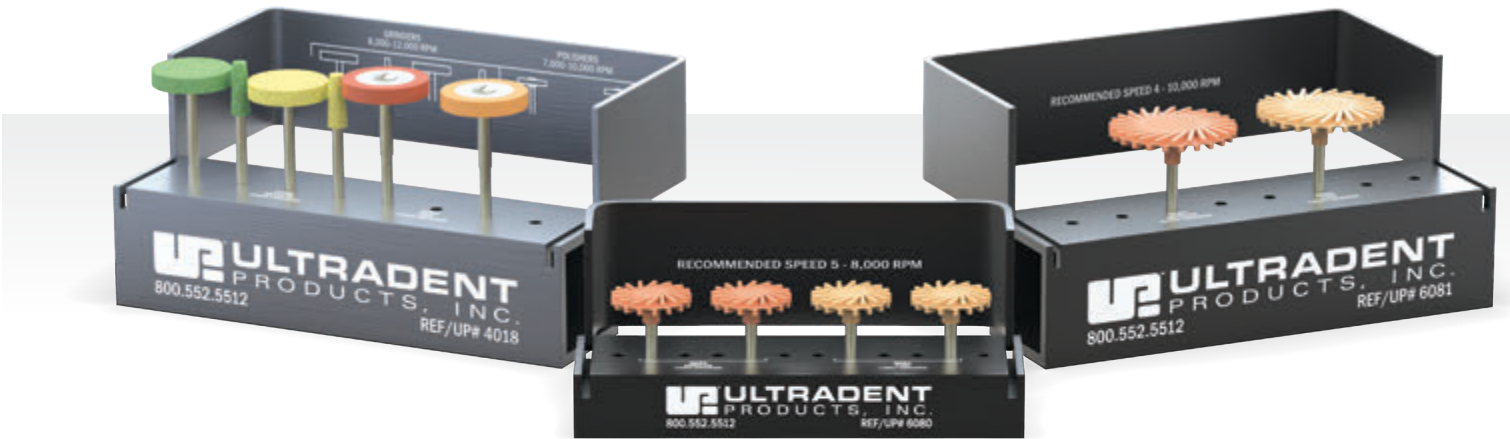


## DON'T SHORTCUT YOUR FINISH

As a dentist, you spend time and effort crafting quality, natural-looking restorations. But if an incorrect or inferior finishing adjuster or polisher is used, the integrity of your restoration can be compromised.

Because of the extreme hardness of newer ceramic materials, effectively adjusting and polishing ceramic restorations can be a major challenge. If a restoration is not adjusted and polished correctly, it can result in microfractures and abrasion of the softer natural dentition of adjacent teeth.

Jiffy Universal Ceramic Kits offer you a universal solution with an optimal number of components, so you can efficiently adjust and polish any ceramic material to achieve a truly smooth surface—with the highest quality, most natural-looking finish possible.

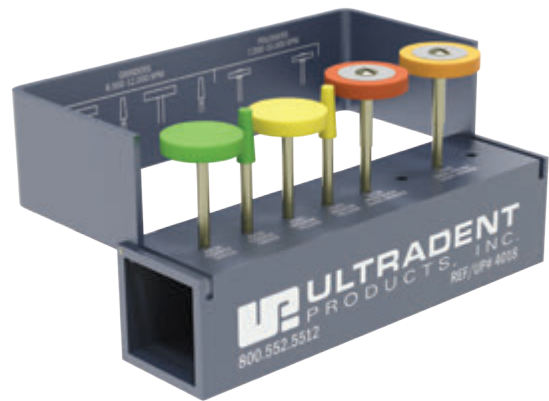

## NATURAL. UNIVERSAL. BEAUTIFUL.

The Jiffy Natural Universal Ceramic Polishing system consists of malleable spiral-shaped wheels that are designed to easily conform to tooth anatomy. Their pliable finger-like extensions easily reach where cups and points can't, and soften super high-gloss finishes for a natural enamel-like result.

They are designed to be used in conjunction with the Jiffy™ Universal Ceramic Adjusting and Polishing System to easily achieve a natural, high-quality finish on all ceramic restorations. The diamond-impregnated wheels are available in extraoral and intraoral wheels in both medium and fine grits. The efficient two-step process allows you to easily achieve the most natural finish on any ceramic material including zirconia, lithium disilicate, and porcelain.

- Naturally adapts to any tooth surface, including occlusal anatomy
- Specially formulated Ultradent diamond grit allows for efficient polishing on any ceramic material, including zirconia
- Designed to temper super high-gloss finishes to achieve a more natural look
- Can be used to restore implant cases
- Designed for use with the Jiffy Universal finishing system
- Efficient two-step process
- Available with or without autoclavable aluminum blocks\*

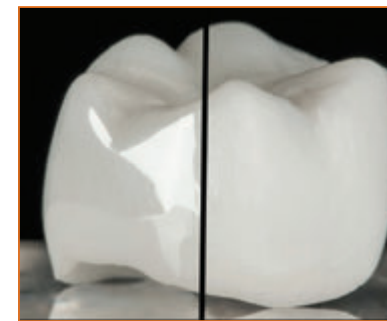


Photo courtesy of Dr. Richard Tuttle, DDS

### What Doctors Are Saying


"These are by far the **quickest** and **best** polishing kits I have ever used. The polish on zirconia is incredible. I have never seen a better, more glass-like appearance on zirconia than I have with the Jiffy Universal System."

—Dr. Eric Hollander DDS, Austin, TX

"These polishing kits give a very nice sheen and are better than any others I've used."

—Dr. William J. Dunbar DDS, Glencloe, MN

### Surface Roughness (uin)<sup>1</sup>

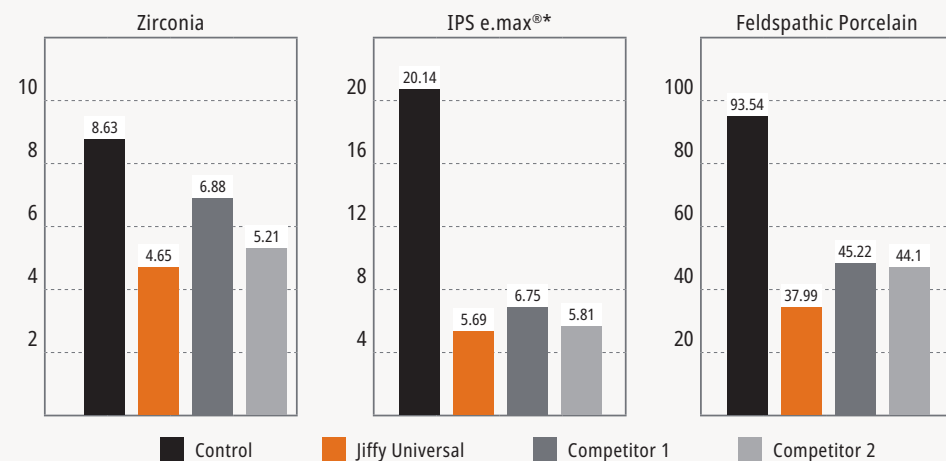

## FEATURES & BENEFITS

- Universal application eliminates the need for multiple adjusters and polishers, saving you time and money
- Specially formulated Ultradent diamond grit provides optimal smoothness and outstanding polishing results, while still being gentle on processed materials
- Multi-grit diamond particles allow for effective adjustment of ceramics for a truly smooth and high-gloss finish
- Maximum diamond retention ensures a long service life
- Autoclavable aluminum block extends the life of the system



### 6080-1 - Jiffy Natural Universal Intraoral Polishing Kit

2 x Jiffy RA Medium Natural Universal 14 mm wheels  
2 x Jiffy RA Fine Natural Universal 14 mm wheels

### 6081-1 - Jiffy Natural Universal Extraoral Polishing Kit

1 x Jiffy HP Medium Natural Universal 26 mm wheel  
1 x Jiffy HP Fine Natural Universal 26 mm wheel

### 6082-1 - Jiffy Natural RA Variety 14 mm Polishing Wheel 6pk

### 6083-1 - Jiffy Natural RA Medium 14 mm Polishing Wheel 3pk 6084-1 - Jiffy Natural RA Fine 14 mm Polishing Wheel 3pk

### 6085-1 - Jiffy Natural HP Medium Polishing 26 mm Wheel 1pk 6086-1 - Jiffy Natural HP Fine Polishing 26 mm Wheel 1pk

<sup>1</sup>. Data on file.  
\*Not a registered trademark of Ultradent Products, Inc.

\*Ultradent recommends the use of a block when autoclaving to prevent warping and deformation