

VALO™ X

BROADBAND LED CURING LIGHT



THE CURING LIGHT
REIMAGINED

NEW!

COMPLETELY REDESIGNED

- Rebuilt from the ground up to create the most innovative curing light available

HIGH-GRADE AEROSPACE ALUMINUM

- Unibody design is exceptionally durable and allows for excellent thermal management

INCREASED LENS SIZE

- 12.5 mm lens covers any tooth while maintaining accessibility and patient comfort

SIMPLIFIED INTERFACE

- Curing and diagnostic modes are indicated, operated, and activated with top and bottom buttons or by accelerometer function

MULTI-CONFIGURATION

- Can be used in a corded or cordless configuration (battery and cord adapter included in kit)

5-YEAR WARRANTY

- Includes a 5-year Manufacturer Warranty

ACCELEROMETER FUNCTION

- Allows you to quickly and easily move through curing and diagnostic modes

12 LED CHIPSET

- Provides high-intensity, broadband light for excellent beam uniformity, curing depth, and beam collimation

CURING MODES: Standard Power Mode, Xtra Power Mode | **DIAGNOSTIC LIGHT MODES:** White Light Diagnostic Aid Mode, Black Light Diagnostic Aid Mode

CYCLING BETWEEN MODES




CURING MODES:
Move the VALO X light forward in a drum tap motion to access and cycle through curing modes.



DIAGNOSTIC LIGHT MODES:
Move the VALO X light in a drum tap motion to the side to access and cycle through diagnostic light modes.

INCLUDED LENS ACCESSORIES


TWO CURING LENSES

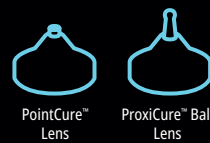

PointCure™ Lens

ProxiCure™ Ball Lens

THREE DIAGNOSTIC AID LENSES

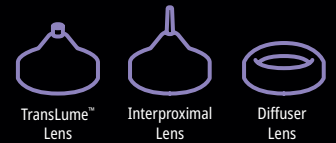

TransLume™ Lens

Interproximal Lens

Diffuser Lens

5973 - VALO X Kit

- 1 x VALO X LED curing light
- 5 x Accessory lenses
- 2 x Rechargeable batteries
- 1 x Battery charger
- 1 x Power supply (for battery charger or cord adapter)
- 1 x Cord adapter
- 1 x Handpiece bracket holder
- 1 x Blue light blocking light shield
- 1 x Sample pack of barrier sleeves

Accessories

- 4951** VALO X Battery Charger
- 4952** Power Supply with Universal Plugs
- 5189** VALO X Assembled Cord
- 5437** VALO X Batteries 2pk
- 4665** VALO X Barrier Sleeves 100pk



SCAN
TO VIEW
PRICING



© 2024 Ultradent Products, Inc. All rights reserved.
1011267-PKAR01 103023

800.552.5512 | ULTRADENT.COM