






# BETTER ACCESSORIES

In most practices, frequently changing tips is part of the routine. The Ultrawave Color-Coding System and wide variety of tips make instrument management simple and efficient.



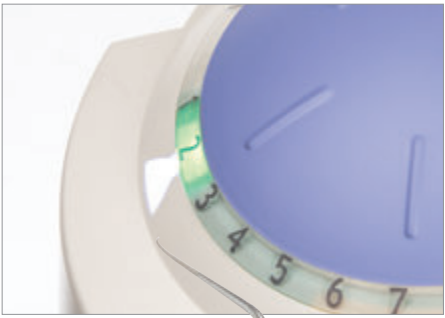



**GREEN**  
Low power, low amplitude.  
Delicate treatments:  
**Periodontics**

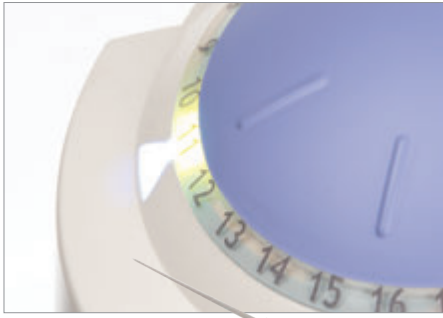





**YELLOW**  
Medium power, medium amplitude.  
Precision treatments:  
**Endodontics & Apical Surgery**

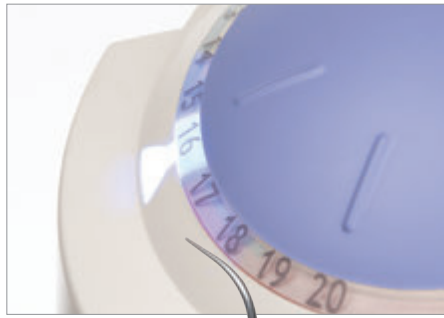

**BLUE**  
High power, high amplitude.  
Routine prophylaxis:  
**Scaling & Restoration**

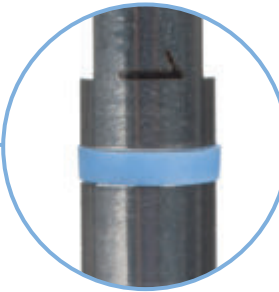

## COLOR-CODING SYSTEM FOR EASY INSTRUMENT MANAGEMENT

The Ultrawave™ and Ultrawave™ XS piezoelectric units offer exceptional ease of use through the Color-Coding System. Matching the right tip to the right power is simple and reliable. The colors green, blue, and yellow allow tips to be used with the most efficient corresponding power setting on the control dial of the devices.

SCALING  
ENDO  
PERIO  
RESTORATIVE

You want an ultrasonic piezo device that saves you time and makes treatment more comfortable for you and your patients. The Ultrawave and Ultrawave XS LED piezo ultrasonic devices do just that. In most practices, frequently changing tips is part of the routine. The Ultrawave Color-Coding System and wide variety of tips makes instrument management simple and stress-free.

# THE ULTRAWAVE TIP ADVANTAGE

- Corrosion resistant: steel alloy is resistant to decontamination products as well as irrigation solutions
- Variety of applications: tips are available for scaling, endodontics, periodontal treatment, and restorative procedures
- Tensile strength: each tip is designed to meet the specific requirements of the procedure
- Sterilization: tips are suitable for autoclavable sterilization




**SCALING**  
Supragingival and subgingival tips are ideal for standard scaling, along with cases where high power is needed for tenacious calculus removal.



**PERIODONTAL**  
Perio tips allow for a full-mouth periodontal debridement. Deposits are removed by methodically crisscrossing all root surfaces using low lateral pressure combined with very low amplitude. This provides tactile sensitivity unequalled by any other ultrasonic tip for more comfortable treatment.



**ENDODONTICS & APICAL SURGERY**  
These tips are specially designed for all the challenging applications during a root canal treatment, including removal of calcification in the pulp chamber, removal of separated instruments, and retreatment. The innovative tip design improves passive ultrasonic irrigation, facilitating superior performance in root canal cleanings.

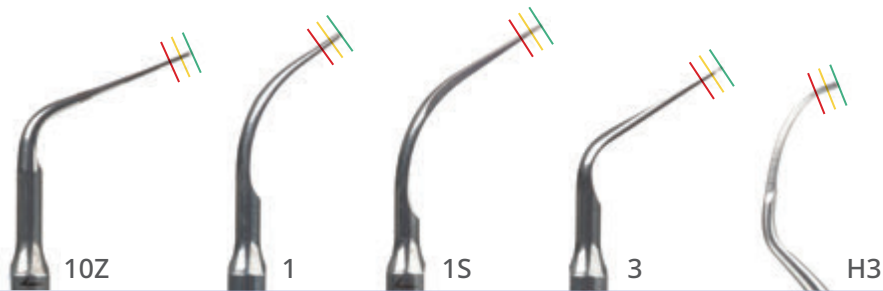

## Ultrawave™ Acteon®\* Satelec®\* Tips

	For supra- and subgingival scaling. Designed for greater power. <b>F00245 - Ultrawave 1 Slim Scaling Tip 1pk</b>		For removal of biofilm and soft deposit. Tip is etched in 3 mm increments. <b>F00254 - Ultrawave 10Z Scaling Tip 1pk</b>
	Recommended for treating simple cases and gross supragingival scaling. <b>F00246 - Ultrawave 1 Scaling Tip 1pk</b>		For removing discolorations and stains. <b>F00248 - Ultrawave 3 Scaling Tip 1pk</b>
	For treating implants, veneers, and porcelain crowns. <b>F00702 - Ultrawave PH1 Tip 4pk</b>		Specially designed for irrigation. <b>F43760 - Ultrawave NiTiSonic Tips 4pk</b>
	Recommended for examining and maintenance of moderate to deep pockets or furcations. <b>F01001 - Ultrawave TK1-1S Tip 1pk</b>		Apical surgery tip. <b>F00065 - Ultrawave AS3 Tip 1pk</b>
	Most popular tip for treatment of anterior teeth. The guiding edge is used parallel to the pocket. <b>F00369 - Ultrawave H3 Tip 1pk</b>		Removes the cement interface between post and canal wall. <b>F88011 - Ultrawave ET20 Tip 1pk</b>
	Left-angled micro tip recommended for maintenance of moderate to deep pockets and furcations. <b>F02161 - Ultrawave TK2-1R Tip 1pk</b>		Tapered smooth tip for retreatment and instrument fragment removal. <b>F88018 - Ultrawave ET25 Tip 1pk</b>
	Right-angled micro tip recommended for maintenance of moderate to deep pockets or furcations. <b>F02162 - Ultrawave TK2-1L Tip 1pk</b>		Diamond-coated tip locates and opens calcified canals. <b>F88020 - Ultrawave ETBD Tip 1pk</b>

## Ultrawave™ Acteon®\* Satelec®\* Wrenches

	<b>F00446 - Flat Wrench</b> <b>F81302 - Blue Wrench</b>
--	------------------------------------------------------------

\* Acteon® and Satelec® are trademarks of Satelec, a company of Acteon Group.