



ExperTemp[®]

Temporary Crown & Bridge Material

Even a temporary restoration should maintain both the function and esthetics of natural dentition. With ExperTemp temporary crown and bridge material you can create provisionals that will closely resemble the final restoration and keep your patients smiling between appointments.

ExperTemp material is a bis-acryl composite used to fabricate temporary crowns, bridges, inlays, and onlays, as well as long-term temporaries. The 10:1 self-cured chemistry is easy to work with and provides exceptional strength, flexibility, and abrasion resistance. Superior performance combined with a natural esthetic make ExperTemp material the ideal choice for temporization.



Are unattractive temporaries making your patients *afraid to smile*?

ExperTemp[®]

Temporary Crown & Bridge Material



ULTRADENT
PRODUCTS, INC.

800.552.5512 ULTRADENT.COM



- Trims easily and polishes beautifully
- Exhibits low oxygen inhibition at polymerization
- Possesses high abrasion resistance and flexibility
- Comes in multiple shades: A1, A2, A3, A3.5, B1, and Bleach White
- Can be repaired with a packable composite, a flowable composite, or additional ExperTemp material
- Fluoresces similarly to enamel



1. ExperTemp temporary crown and bridge with esthetic translucency



2. Prepped teeth



3. Etch (2-secs) removes film from hemostatic agent



4. ClearTemp® LC temporary cement applied



5. Blade used to trim temporary crown and bridge without altering margins




6. ExperTemp achieves esthetic blend with natural teeth


Technical Data¹

Compressive strength	447 MPa
Flexural strength	127 MPa
Hardness (Durometer)	89 points
Volume Shrinkage	5.1%
Working time	40 sec.


Edge Chip

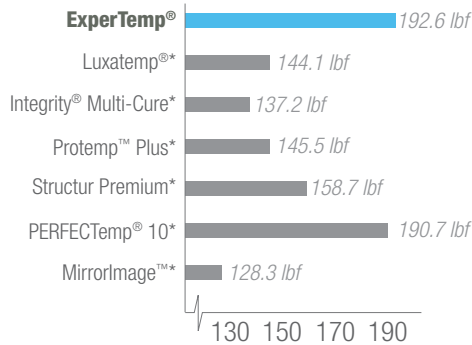

Flex Modulus

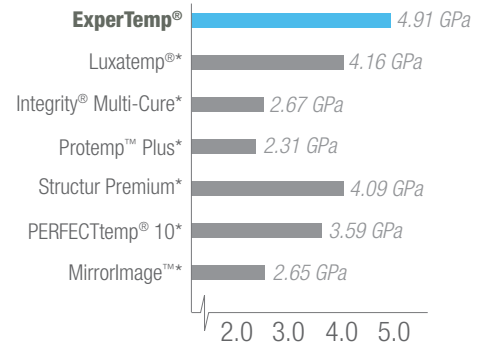

*Trademark of a company other than Ultradent. 1. Data on file

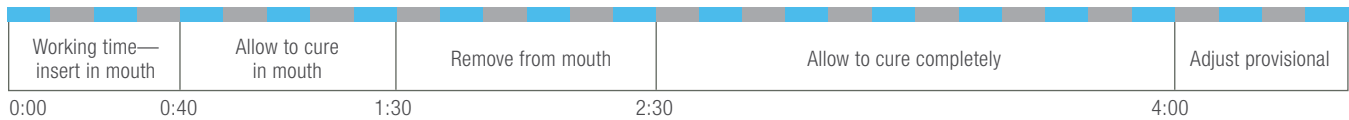

- 6340 - ExperTemp A2 Kit
- 6341 - ExperTemp A1 Kit
- 6342 - ExperTemp A3.5 Kit
- 6343 - ExperTemp B1 Kit
- 6344 - ExperTemp Bleach White Kit
- 6347 - ExperTemp A3 Kit

- 6345 - ExperTemp Dispensing Gun
1 x 10:1/4:1 dispensing gun

- 6346 - ExperTemp Mixing Tip Refill 45pk
45 x Mixing tips



All kits include one 50 ml cartridge and 15 mixing tips