

VISCOSTAT™ HEMOSTATIC




3093 - ViscoStat Dento-Infusor Syringe Kit
4 x 1.2 ml syringes
20 x Metal Dento-Infusor tips



647 - ViscoStat Dento-Infusor IndiSpense Syringe Kit
1 x 30 ml IndiSpense syringe
20 x Metal Dento-Infusor tips
20 x 1.2 ml empty syringes



3070 - ViscoStat Syringe 20pk
1.2 ml syringes

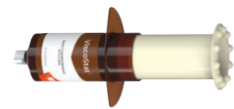


645 - ViscoStat IndiSpense Syringe 1pk
30 ml syringe




1278 - Hemostatic Delivery Syringe 20pk
3096 - Hemostatic Delivery Syringe 100pk
1.2 ml empty syringes

ASTRINGEDENT™ X HEMOSTATIC

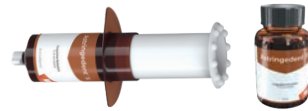

112 - Astringedent X Bottle 1pk
690 - Astringedent X IndiSpense Syringe 1pk
30 ml

ASTRINGEDENT™ HEMOSTATIC

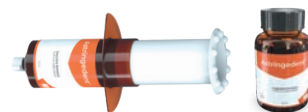

111 - Astringedent Bottle 1pk
686 - Astringedent IndiSpense Syringe 1pk
30 ml


METAL DENTO-INFUSOR™ TIP WITH COMFORT HUB™ DESIGN

LOK-TITE + COMFORT HUB	100pk	500pk
19 ga Metal Dento-Infusor	4955	4956



BLUE MINI™ DENTO-INFUSOR™ TIP

LOK-TITE	mm	100pk	500pk
Blue Mini Dento-Infusor	1.20	1086	1440

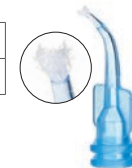

! Note: Do not use epinephrine preparations with ferric sulfate products (ViscoStat, Astringedent, Astringedent X hemostatics), as blue/black precipitate will occur.



SCAN TO VIEW PRICING

FOR PROFOUND HEMOSTASIS



Tissue Management

Unparalleled tissue management starts with rapid, profound hemostasis. For more than 40 years, dentists have trusted the immediate hemostatic power, detailed margins, and elimination of surface bleeding and sulcular fluid provided by Ultradent's tissue management products.

Technique Guide



1 Bleeding at preparation.



2 Astringedent X hemostatic is scrubbed firmly against bleeding tissue.



3 When bleeding stops, give a firm air/water spray to test for hemostasis and remove coagulum.



4 Profound hemostasis.



ASTRINGEDENT™ X HEMOSTATIC

12.7% Iron Solution

The preferred hemostatic for most case situations, not just difficult ones

- Ultradent's fastest and most powerful hemostatic*
- Stops minor to severe bleeding



ASTRINGEDENT™ "THE CLASSIC" HEMOSTATIC

15.5% Ferric Sulfate

Known as the "classic" hemostatic agent—Ultradent's first product

- Eliminates sulcular fluid contamination for optimal bonding



VISCOSTAT™ HEMOSTATIC

20% Ferric Sulfate

Viscous formula allows for superior control


- Gentle on tissues
- Creates rapid, profound hemostasis
- Viscous solution won't run, drip, or migrate to other teeth

*Data on file.

TEST YOUR SCRUBBING SKILLS

Are you using adequate pressure to achieve rapid, profound hemostasis?

Directions: Use your Dento-Infusor™ tip to scrub the scratch-off area.

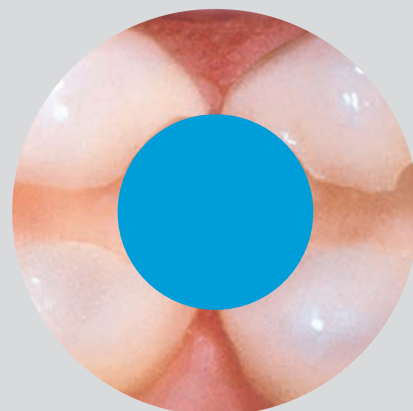

- Metal Dento-Infusor™ tip with Comfort Hub™ design has a metal, bent cannula with a flared brush padded end. The larger ergonomic wings provide a secure, comfortable grasp of the tip.
- Blue Mini™ Dento-Infusor™ tip has a softer padded end
- Designed to work with all our hemostatic agents
- Designed to be a tiny scrub brush, not only an applicator
- Facilitates infusion of hemostatics and scrubbing tissue clean of coagulum