



VALO X Colors Kit –

- 1 x VALO X LED curing light
- 1 x Diffuser lens
- 2 x Rechargeable batteries
- 1 x Battery charger
- 1 x Blue-light-blocking shield
- 1 x Sample pack of barrier sleeves
- 1 x Power supply (for battery charger or cord adapter*)

- **6496** Red
- **6497** Pink
- **6498** Gray
- **6499** Green
- **6519** Blue
- **6538** Purple
- **6539** Black

5973 - VALO X Kit –

- 1 x VALO X LED curing light
- 5 x Accessory lenses
- 2 x Rechargeable batteries
- 1 x Battery charger
- 1 x Power supply (for battery charger or cord adapter)
- 1 x Cord adapter
- 1 x Handpiece bracket holder
- 1 x Blue-light-blocking shield
- 1 x Sample pack of barrier sleeves



6383 - VALO XG Kit –

- 1 x VALO XG LED curing light
- 2 x Rechargeable batteries
- 1 x Battery charger
- 1 x Blue-light-blocking shield
- 1 x Sample pack of barrier sleeves
- 1 x Power supply (for battery charger or cord adapter*)

6470 - VALO XS Kit –

- 1 x VALO XS LED curing light
- 2 x Rechargeable batteries
- 1 x Battery charger
- 1 x Blue-light-blocking shield
- 1 x Sample pack of barrier sleeves
- 1 x Power supply (for battery charger or cord adapter*)

Accessories


- 4951** Battery Charging Unit
- 4952** Power Supply
- 5189** Assembled Cord
- 5437** Rechargeable Batteries *2pk*
- 4665** Barrier Sleeves *100pk*
- 5898** Barrier Sleeves *600pk*



5839 Curing Light Ease-In-Shields Kit

Lens Accessories


TWO CURING LENSES

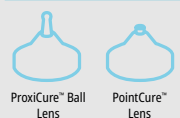

5889 (X)
4081 (XG)
5936 (XS)

5890 (X)
4082 (XG)
5934 (XS)

FOUR DIAGNOSTIC AID LENSES

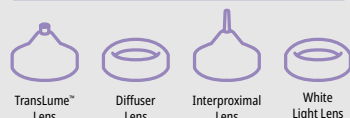

5891 (X)
4084 (XG)
5937 (XS)

5892


5893 (X)
4658 (XG)
4629 (XS)

4628

WELCOME TO THE FAMILY

The VALO X curing light family is expanding.

Built on the same foundation of power, accessibility, and durability, the new VALO XG and VALO XS curing lights join the VALO X light to offer multiple models, so you can choose the curing light that best fits your office.

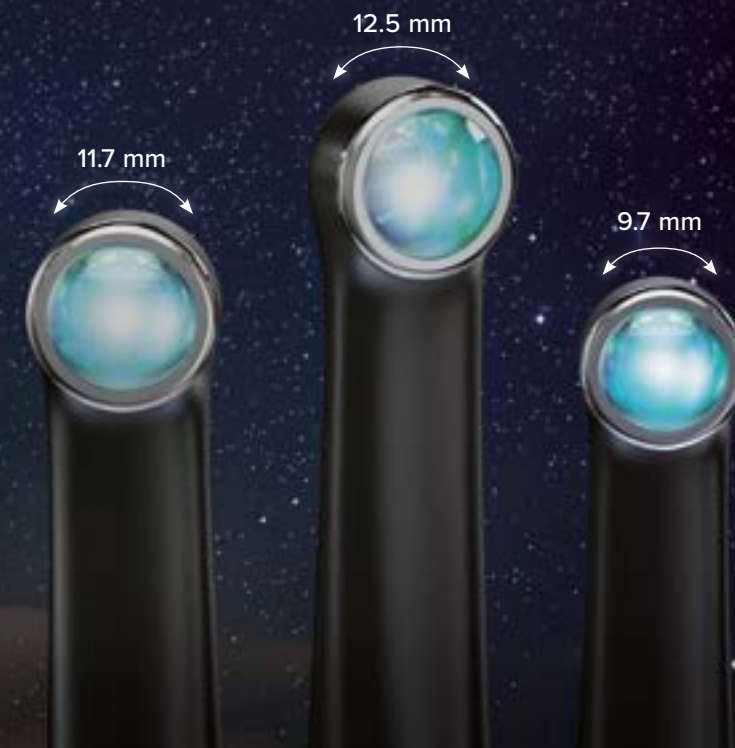

*Cord adapter not included



SCAN
TO VIEW
PRICING

DURABILITY YOU LOVE

A curing light is an essential part of your office. It needs to stand up to day-to-day use. If your curing light is out of commission, it can cost you a lot of time, money, and hassle. With the VALO X curing light family, you won't have to worry, because each model can handle it.

All VALO X curing lights are created from a solid bar of high-grade aerospace aluminum, making each model exceptionally durable. The unibody construction creates a sealed, strong light that eliminates weak connection points and allows for excellent heat dissipation.



The VALO X curing light family is made to last, and we stand behind it. Don't worry about having to replace your curing light every couple of years; every VALO X light comes with a 5-year manufacturer's warranty so you can be confident in the value of your purchase.

ACCESSIBILITY YOU WANT

VALO X, VALO XG, and VALO XS lights have a low-profile design to allow for easy access to posterior restorations without sacrificing patient comfort. The slim head allows the light to be placed directly over the curing site, no matter where it is in the mouth, ensuring light can reach all aspects of the preparation. The VALO X light family's larger lenses give a larger curing surface area, so you can get the right amount of light in the right place.



BROADBAND LEDS

- Optimally collimated multi-wavelength beam delivers consistent, uniform power that can cure all light-curable dental materials

ACCELEROMETER FUNCTION

- Allows you to quickly and easily move through curing and diagnostic aid modes

SIMPLIFIED INTERFACE

- Curing and diagnostic aid modes are indicated, operated, and activated with top and bottom buttons or by Accelerometer Function

MULTI-CONFIGURATION

- Can be used in a corded or cordless configuration

A COMPLETE CONSTELLATION OF CURING LIGHTS

VALO X

VALO XG

VALO XS



LENS SIZE	12.5 MM	11.7 MM	9.7 MM
SURFACE AREA	123 MM ²	108 MM ²	74 MM ²
STANDARD CURING MODE	1100 MW/CM ² (10 SECONDS)	1000 MW/CM ² (10 SECONDS)	1000 MW/CM ² (10 SECONDS)
XTRA POWER CURING MODE	2200 MW/CM ² (5 SECONDS)	2000 MW/CM ² (5 SECONDS)	2000 MW/CM ² (5 SECONDS)
ORTHO QUADRANT CURING MODE	—	—	2000 MW/CM ² (5 X 3 SECONDS)
BUILT-IN DIAGNOSTIC AID MODES	BLACK LIGHT (1 MIN.) WHITE LIGHT (1 MIN.)	BLACK LIGHT (1 MIN.)	BLACK LIGHT (1 MIN.)
ACCESSORY LENS COMPATIBILITY	✓	✓	✓
DURABLE UNIBODY DESIGN	✓	✓	✓
BROADBAND LEDS	✓	✓	✓
CORDED/CORDLESS OPERATION	✓	✓	✓
ACCELEROMETER	✓	✓	✓
5-YEAR WARRANTY	✓	✓	✓
12 MM MOLAR COVERAGE			
BEAM PROFILE			