



Our NaviTip™ tips are treated to create flexible tips with a smooth, round end—ensuring they will easily snake through canals to deliver irrigants precisely where they are needed. This treatment can create a golden hue you may see on some of our NaviTip tips.

NaviTip™ Tips 29 ga with Single Sideport

- Designed to direct the flow of chemistry through the sideport of the tip before flowing down into the area of the apex, thus reducing the risk of product extrusion
- Flexible, stainless steel cannula easily navigates curved canals
- Bendable tip 0.33 mm**

29 ga delivers paste materials: MTAFlow™, MTApeX™, and UltraCal™ XS.

Note: UltraCal™ XS calcium hydroxide paste should only be used with NaviTip 29 ga Single Sideport tips.

| | Tip length | 20pk |
|-------------------------|------------|------|
| 29 ga - 0.0129" NaviTip | 27 mm | 4989 |
| 29 ga - 0.0129" NaviTip | 25 mm | 4990 |
| 29 ga - 0.0129" NaviTip | 21 mm | 4991 |
| 29 ga - 0.0129" NaviTip | 17 mm | 4992 |
| 29 ga-29 ga NaviTips | 27-17 mm | 5143 |

NaviTip™ Tips 31 ga with Double Sideport Irrigator

- Double sideports deliver irrigants safely, minimizing the possibility of chemicals being expressed past the apex
- The world's smallest cannula navigates the most intricate canal spaces

Designed for: ChlorCid™, ChlorCid™ Surf, Ultradent™ EDTA 18% Solution, Consepis™, and Citric Acid 20%.

| | |
|-------------------|-------------------|
| 5121 - 21 mm 20pk | 5123 - 27 mm 20pk |
| 5122 - 21 mm 50pk | 5124 - 27 mm 50pk |

NaviTip™ Tips 29 ga and 30 ga

- Provide controlled delivery to the apex
- Flexible, stainless steel cannula easily navigates curved canals
- Available in 17, 21, 25, and 27 mm lengths and 29 or 30 ga cannulae

29 ga delivers paste materials.

30 ga delivers solutions and gels.

| | 29 ga - 20pk | 30 ga - 20pk | 29 ga - 50pk | 30 ga - 50pk |
|-------------------|--------------|--------------|--------------|--------------|
| 17 mm - White | 5112 | 1249 | 1378 | 1421 |
| 21 mm - Yellow | 5113 | 1349 | 1374 | 1422 |
| 25 mm - Blue | 5114 | 1250 | 1376 | 1423 |
| 27 mm - Green | 5115 | 1354 | 1377 | 1424 |
| Assorted 17-27 mm | 5116 | 1351 | 1379 | 3319 |

LOK-TITE™

new



LOK-TITE™

NaviTip™ FX™ Tips 30 ga


- One-of-a-kind brush cleans, scrubs, and irrigates simultaneously
- Rigid 30 ga cannula

Designed for: Ultradent™ Citric Acid 20% Solution.

Listed as an "EXCELLENT" product by a prominent independent research institute.¹

| |
|-------------------|
| 1452 - 17 mm 20pk |
| 1454 - 25 mm 20pk |

1. Clinical Research Associates Newsletter. Volume 29, Issue 1, January 2005.

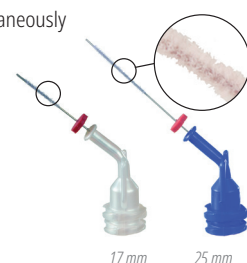

Endo-Eze™ Irrigator Tip 27 ga

- Provides ideal reach reducing risk of expressing chemicals past the apex
- Comes with a flexible, blunt cannula with a unique, anti-obturator end
- Non-sterile

Designed for: Ultradent™ 5 ml syringe.

| |
|------------------|
| 207 - 25 mm 20pk |
|------------------|



Endo-Eze™ Tips

- Work for a variety of endodontic procedures
- Come with flexible, strong cannulae in a variety of sizes
- Bend easily

Designed for: luting materials and air/water delivery.

Use with: TriAway™ Adapter, PermaFlo™ DC (20 ga), and other Ultradent syringes.

| Gauge | Diameter | 20pk | 100pk |
|-------|-----------------|------|-------|
| 18 ga | 0.049"; 1.25 mm | 345 | 1428 |
| 19 ga | 0.042"; 1.06 mm | 346 | 1429 |
| 20 ga | 0.035"; 0.90 mm | 347 | 1430 |
| 22 ga | 0.028"; 0.70 mm | 348 | 1431 |



Capillary Tips

- Evacuates canals and substantially minimizes use of paper points
- Narrow, flexible taper accesses curved canals
- Do not use to deliver irrigating materials or endodontic chemistries**

Designed for: Ultradent™ Luer Vacuum Adapter for moisture removal from endodontic canals.

| Diameter | 20pk | 50pk |
|-----------------|------|------|
| 0.014"; 0.36 mm | 341 | 3099 |
| 0.019"; 0.48 mm | 186 | 1425 |

LOK-TITE™



Ultradent Tips

DESIGNED TO DELIVER



SCAN TO VIEW PRICING



The chemistries you use are different. Some are chemically activated—needing to be mixed immediately before delivery. Others have varying viscosities. Some work in pits and fissures, some inside canals, and some on enamel surfaces. Each chemistry you use is designed for a specific purpose. Shouldn't the same be true for your tips?

Ultradent makes tips **designed to deliver** each chemistry we create. Whether you're delivering a flowable composite, a viscous gel, or thick impression material we make the perfect tip for the job. And since our tips are engineered on-site, we test each design to ensure it works perfectly with the chemistry it's intended for.

Check out our tips with **LOK-TITE™** and **COMFORT HUB™**

Luer Lock tips with Lok-Tite feature double threads that lock the tip into place for increased security and wings for easy attachment and removal. Tips with the Comfort Hub feature include larger ergonomic wings that provide a secure, comfortable grasp of the tip.

Black Mini™ Tip

- Dispenses large volumes
- Opaque plastic preserves the flow of light-cured materials

Designed for: Ultra-Blend™ plus, Ultradent™ LC Block-Out Resin, PermaFlo™, TriAway Adapter, UltraTemp™, Opalescence™ Boost™, Ultradent™ Diamond Polish Mint, OpalDam™, Opalescence™ Endo, and OraSeal™ Caulking.

| |
|--------------|
| 514 - 100pk |
| 1433 - 500pk |

LOK-TITE™



Black Mini™ Brush Tip

- Precise, controlled delivery of aqueous materials
- Tight, adjustable brush fibers minimize bubbles
- Unique to Ultradent

Designed for: Consepsis™, Peak™ SE Primer, Peak™-ZM, Seek™/Sable™ Seek™, Ultradent™ Silane, and Ultradent™ Universal Dentin Sealant.

| |
|--------------|
| 1169 - 100pk |
| 1432 - 500pk |

LOK-TITE™



Black Micro™ FX™ Tip

- Accommodates various viscosities
- Flocked tip fans out to spread materials in a thin, uniform layer
- 22 ga tip

Designed for: Composite Wetting Resin, PQ1, PrimaDry™ and PermaSeal™.

| |
|--------------|
| 1357 - 100pk |
| 1434 - 500pk |

LOK-TITE™



Black Micro™ Tip · Blue Micro™ Tip

- Pinpoint precision
- Narrow cannulae accurately deliver materials
- Black Micro tip: 22 ga, Blue Micro tip: 25 ga

Black Micro tip designed for: Ultra-Blend™ plus and EtchArrest™. Blue Micro tip designed for: Ultra-Etch™, PrimaDry™, and Opal-Etch™.

| | |
|--------------------|-------------------|
| 1085 - Black 100pk | 127 - Blue 100pk |
| 1435 - Black 500pk | 1436 - Blue 500pk |



Intraoral Impression Tip

- Allows precise placement
- Attaches to Impression Mixing tip

Designed for: Thermo Clone™ VPS and Chromaclone™ VPS.

| |
|-------------|
| 2906 - 50pk |
|-------------|



Impression Mixing Tips

- Enables direct delivery of impression materials
- Automixing, disposable, and color coded

Designed for: Thermo Clone VPS and Chromaclone VPS.

| |
|--------------------|
| 2902 - Yellow 50pk |
| 2903 - Pink 50pk |
| 2904 - Green 50pk |



SoftEZ™ Tip

- Tip fibers provide visible, controlled delivery
- Brush fibers facilitate smooth application

Designed for: Enamelast™.

| |
|--------------|
| 4712 - 50pk |
| 4711 - 100pk |

LOK-TITE™



White Mac™ Tip · White Mini™ Tip

- Dispense large volumes
- All-plastic delivery tips

Designed for: OraSeal™ Caulking, OraSeal™ Putty, Consepsis™ Scrub, and Opalustre™.

| |
|------------------------|
| 1361 - White Mac 100pk |
| 1480 - White Mac 500pk |
| 1247 - White Mini 20pk |



Intraoral Tip

- Allows precise placement
- Attaches to dual-barrel mixing tips

Designed for: PermaFlo™ DC.

| | |
|-------------|------------|
| 5922 - 20pk | \$11.99 ea |
|-------------|------------|



LOK-TITE™

Micro 20 ga Tip

- Large-gauge cannula enables consistent flow
- Standard flowable composite delivery tip

Designed for: MTAFlow™, PermaFlo, PermaFlo™ Purple, OpalDam, and OpalDam™ Green.

| |
|--------------|
| 1252 - 100pk |
| 1437 - 500pk |



COMFORT HUB™ LOK-TITE™

Metal Dento-Infusor™ Tip with Comfort Hub™

- Comfort Hub™ design for secure, comfortable control while attaching or removing
- Places hemostatic agents precisely and effectively removes superficial coagulum
- Blunt, bent 19 ga cannula with padded brush enables gentle pressure on the sulcus
- Ultradent's first tip, the "MDI" remains paramount for successful tissue management

Designed for: Astringedent™, Astringedent™ X, ViscoStat™, ViscoStat™ Clear, and Peak™ Universal Bond.

| |
|--------------|
| 4955 - 100pk |
| 4956 - 500pk |




Blue Mini™ Dento-Infusor™ Tip

- Offers the same tissue management benefits as the Metal Dento-Infusor tip
- Allows controlled flow of drop-sized quantities
- Outer diameter 1.2 mm

Designed for: Hemostatics and Consepsis.

| |
|--------------|
| 1086 - 100pk |
| 1440 - 500pk |

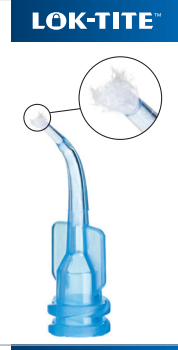

LOK-TITE™

Metal Opal™ Seal Tip

- Soft fibers are ideal for applying a consistent, even, thin layer to etched enamel
- Metal cannula provides ideal flow and ease of expression

Designed for: Opal™ Seal.

| |
|-----------------|
| 500-094 - 20pk |
| 500-095 - 100pk |
| 500-097 - 500pk |



LOK-TITE™

Ultradent™ Mixing Tip

- Mixes and delivers in one action

Designed for: UltraTemp™, UltraTemp™ REZII, PermaFlo™ DC, and EndoREZ™.

| |
|-------------|
| 5920 - 20pk |
|-------------|



SST™ (Surgical Suction Tip)

- Ideal for delicate surgeries
- Large-diameter tip opening

Designed for: Ultradent™ Luer Vacuum Adapter for small periodontic or endodontic procedures and controlled suction of Opalescence Boost.

| |
|-------------|
| 1248 - 20pk |
|-------------|

LOK-TITE™



Inspiral™ Brush Tip

- Delivers viscous or filled materials smoothly via an internal helical channel and ridge
- Tight, adjustable brush fibers minimize bubbles

Designed for: Composite Wetting Resin, Peak Universal Bond, Ultradent™ Porcelain Etch, Ultra-Etch, UltraSeal™ XT plus, and UltraSeal™ XT hydro.

| |
|----------------|
| 123 - 100pk |
| 1033-I - 500pk |

LOK-TITE™



Micro Capillary™ Tips

- Bright colors are easily identified against soft tissues
- The world's smallest molded tips (inside diameter: 0.008")

Designed for: Perio materials, endodontics, and Ultradent™ Luer Vacuum Adapter.

| | |
|-------------------|-------------------|
| 1120 - 5 mm 20pk | 1121 - 10 mm 20pk |
| 1802 - 5 mm 100pk | — |

LOK-TITE™



Opal™ Seal Tip

- Soft fibers are ideal for applying a consistent, even, thin layer to etched enamel
- Specially designed to deliver filled resins with internal helical channel and ridge

Designed for: Opal™ Seal.

| |
|-----------------|
| 500-064 - 20pk |
| 500-034 - 100pk |

LOK-TITE™

