
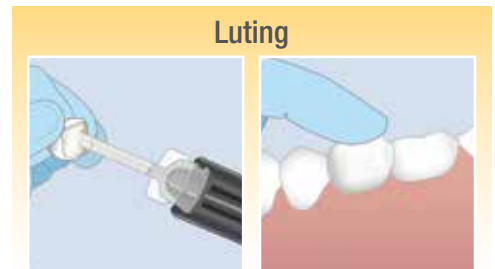
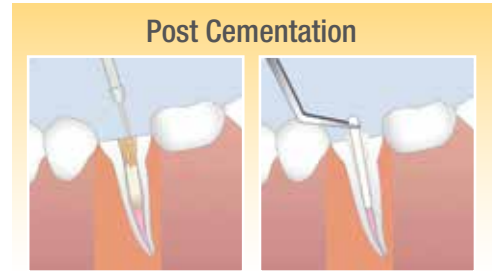


# PermaFlo® DC

COMPOSITE LUTING/RESTORATIVE RESIN

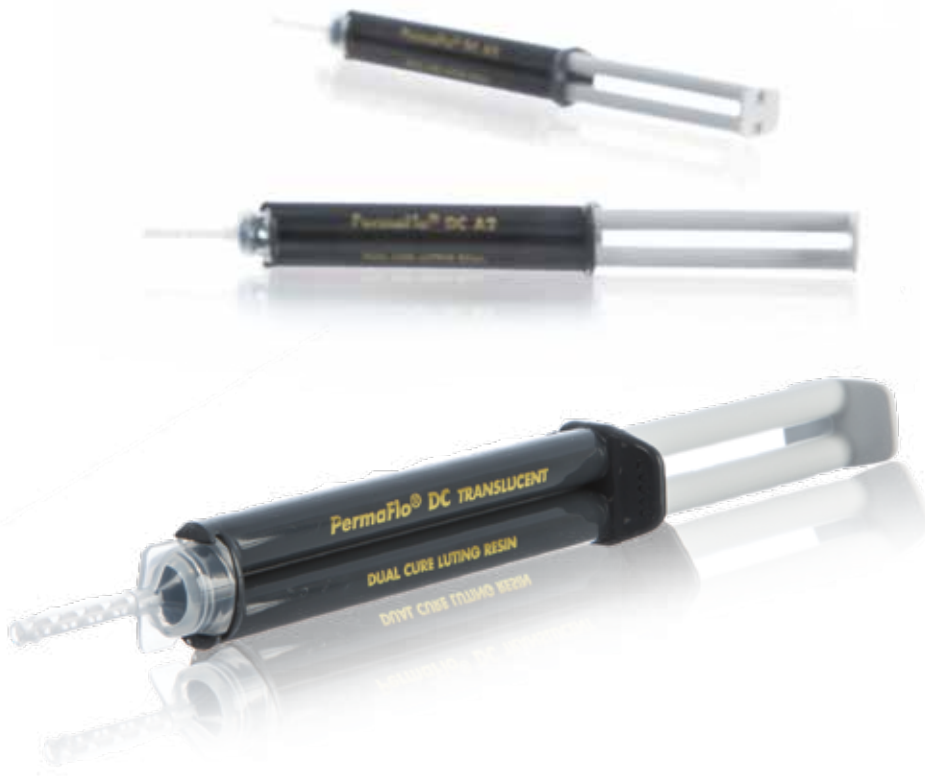

# MULTIPLE OPTIONS

## Failure is NOT one of them

Too often, clinicians see restorations fail sooner than expected. But with PermaFlo DC composite resin, you can give every patient lasting results. This versatile material is compatible with nearly any substrate and can be used for luting, post cementation, or fabricating core buildups. The dual-cure composite resin formula is highly filled, providing the durability and wear resistance for a restoration that will withstand the test of time.

Use PermaFlo DC composite resin to lute inlays, onlays, crowns, and bridges made of any of these compatible substrates:\*

• Porcelain / feldspar-based ceramics	• Leucite-based ceramics
• Lithium disilicate (IPS e.max®)	• Alumina-based ceramics
• Zirconia	• Mica-based ceramics
• Porcelain fused to metal	• Gold
• Metal alloys	• Composite

\*Recommended bonding adhesive is Peak® Universal Bond adhesive



ULTRADENT  
PRODUCTS, INC.

800.552.5512 ULTRADENT.COM



- Thixotropic with low film thickness (9 µm)
- Highly filled – 70% by weight
- Radiopaque with low polymerization shrinkage
- Compatible with total-etch or self-etch bonding systems
- Available in four shades:  
A2, A3.5, Translucent, Opaque White

## Technical Data<sup>1</sup>

Shear Bond Strength to Enamel (Total-Etch)	53.38 MPa
Shear Bond Strength to Dentin (Total-Etch)	62.07 MPa
Flexural Strength	128.5 MPa
Flexural Modulus	9.37 GPa
Compressive Strength	355.91 MPa
Compressive Modulus	4.22 GPa

1. Data on file.



- 1319 PermaFlo DC Bonded Post Kit**
- 2 x 5 ml dual-barrel PermaFlo DC A2 syringes
  - 2 x 1.2 ml Ultra-Etch syringes
  - 2 x 1.0 ml Peak SE syringes
  - 2 x 1.2 ml Peak Universal Bond syringes
  - 5 x Luer Vacuum Adapters
  - 20 x TriAway Adapters
  - 20 x Skini syringes
  - 20 x Intraoral tips
  - 20 x Mixing tips
  - 40 x NaviTip FX 17 mm tips
  - 40 x Endo-Eze 22 ga tips



- 1318 PermaFlo DC Indirect Luting Kit**
- 3 x 5 ml dual-barrel PermaFlo DC syringes  
(one each of A2, A3.5, and Translucent)
  - 2 x 1.0 ml Peak SE syringes
  - 1 x 1.2 ml Consepsis Scrub syringe
  - 2 x 1.2 ml Peak Universal Bond syringes
  - 1 x 1.2 ml Ultra-Etch syringe
  - 20 x Skini syringes
  - 40 x Inspiral Brush tips
  - 20 x Black Mini tips
  - 20 x Mixing tips
  - 20 x Intraoral tips
  - 40 x Black Mini Brush tips
  - 40 x Blue Micro tips
  - 1 x PropGard regular
  - 2 x Tongue guards



- |             |                 |             |                         |
|-------------|-----------------|-------------|-------------------------|
| <b>5912</b> | <b>A2 Kit</b>   | <b>5914</b> | <b>Translucent Kit</b>  |
| <b>5913</b> | <b>A3.5 Kit</b> | <b>5915</b> | <b>Opaque White Kit</b> |
- 1 x 5 ml dual-barrel PermaFlo DC syringe
  - 20 x Mixing tips
  - 20 x Intraoral tips

- 5920 Mixing Tips**  
20 x Mixing tips



- 5922 Intraoral Tips**  
20 x Intraoral tips

