



# THE ORIGINAL KNITTED CORD


with Ultrapak™ CleanCut

---

Easy packing. Excellent absorption. Exceptional retraction.


## Ultrapak™ Knitted Cord

*Available in six versatile sizes without epinephrine and four sizes with epinephrine.*

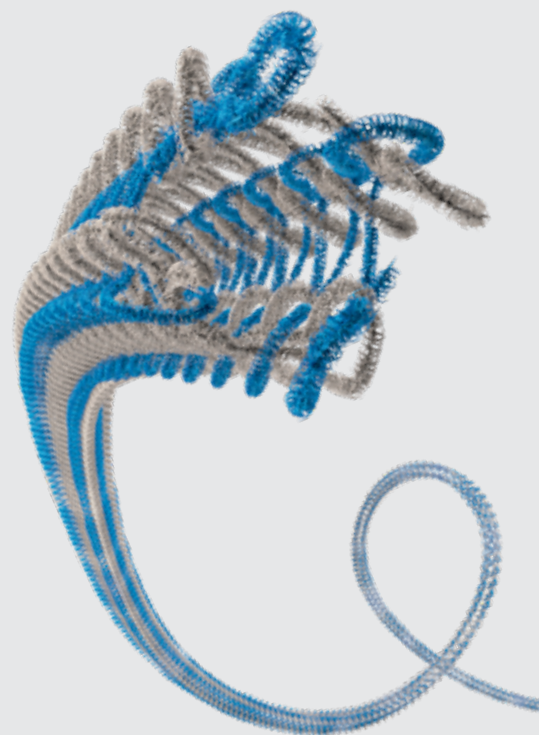


# Ultrapak™ Knitted Cord

When tissue management is key, reach for Ultrapak knitted displacement cord. Its proprietary knitted design is better than braided or twisted cord for **packing** into the sulcus, **absorbing** fluids, and **expanding** for adequate retraction. As an important part of Ultradent's proven tissue management system, Ultrapak cord ensures maximum control with minimum effort.

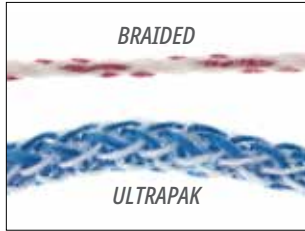

Knitted Ultrapak cord is composed of thousands of tiny, interlocking loops so it compresses and expands easier than other cords. 100% cotton fibers provide high absorption of hemostatic agents and sulcular fluids.

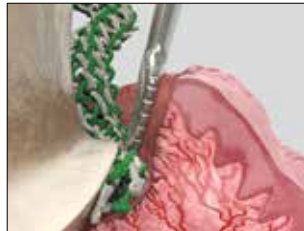

Ultrapak cord packs effortlessly into the sulcus and then exerts a gentle, continuous outward force as knitted loops seek to open. Optimal tissue displacement occurs in 1-3 minutes.



Ultrapak packing cord does not entangle in burs or packers. The bright color facilitates easy location and removal.



Ultrapak CleanCut design features a blade in the cap for efficient cutting. A special dispensing orifice prevents cord from falling into the bottle.

*Ask about our Mix and Match program*

**137 - Ultrapak Cord #000 Regular 1pk**



**132 - Ultrapak Cord #1 Regular 1pk**



**199 - Ultrapak Cord #1 Epinephrine 1pk**



**136 - Ultrapak Cord #00 Regular 1pk**



**133 - Ultrapak Cord #2 Regular 1pk**



**197 - Ultrapak Cord #00 Epinephrine 1pk**

**200 - Ultrapak Cord #2 Epinephrine 1pk**

**131 - Ultrapak Cord #0 Regular 1pk**



**134 - Ultrapak Cord #3 Regular 1pk**



**198 - Ultrapak Cord #0 Epinephrine 1pk**

**130 - Ultrapak Kit**

1 x Each #00, 0, 1, and 2 cord  
1 x Ultrapak organizer



*Each bottle contains 8 ft/244 cm of cord.*

*Ultrapak E is impregnated with racemic epinephrine hydrochloride.*



**Recommended Product** Fischer's Ultrapak Packers  
Specially designed packers ease the packing of Ultrapak knitted cord. Available with thin, serrated or non-serrated edges, they allow practitioners to place the cord in one circular, sliding motion without lifting the instrument from the cord.

