

NaviTip™ Delivery Tips

When it comes to root canals, Ultradent is with you every step of the way.



NaviTip™ Double Sideport Tip

NaviTip™ Tip

NaviTip™ FX™ Tip

Controlled Delivery

Ultradent's NaviTip tips provide the optimal irrigant delivery for every step in a root canal procedure. Unlike traditional irrigation tips, the cannula of each NaviTip tip is slightly rigid through the base and center, but flexible at the tip to allow for easy navigation of curved canals and precise delivery of irrigants short of the apex.

NaviTip™ Tips

LOK-TITE™	20pk		50pk	
	29 ga	30 ga	29 ga	30 ga
17 mm, White	5112	1249	1378	1421
21 mm, Yellow	5113	1349	1374	1422
25 mm, Blue	5114	1250	1376	1423
27 mm, Green	5115	1354	1377	1424
Assorted	5116	1351	1379	3319

NaviTip™ Double Sideport Irrigator Tips

LOK-TITE™	mm	20pk	50pk
31 ga NaviTip	21	5121	5122
31 ga NaviTip	27	5123	5124

NaviTip™ FX™ Tips

LOK-TITE™	mm	20pk
30 ga NaviTip FX	17	1452
30 ga NaviTip FX	25	1454




NaviTip™ 29 ga or 30 ga Tip

- Controlled delivery anywhere in the canal
- Flexible cannula easily navigates curved canals



NaviTip™ 31 ga Double Sideport Irrigator Tip

- The small cannula navigates the most intricate canals
- Double sideports horizontally express irrigant toward the canal walls, not the apex

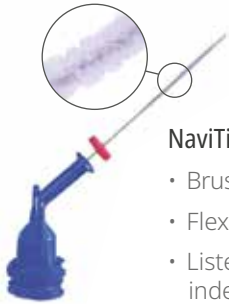


NaviTip™ 30 ga FX™ Tip

- Brush tip scrubs, cleans, and irrigates simultaneously
- Flexible cannula easily navigates curved canals
- Listed as an "EXCELLENT" product by a prominent independent research institute*

NaviTip™ Tips Adapt to Your Style

NaviTip tips deliver low-viscosity materials exactly where they are needed—even in the smallest canals.

The double sideport openings in the NaviTip 31 ga tip allow your irrigant to be safely expressed against the canal walls rather than toward—or even through—the foramen.

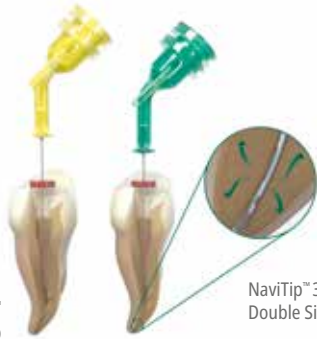

NaviTip™
30 ga Tip

NaviTip™ 31 ga
Double Sideport Tip

*Clinical Research Associates Newsletter, Volume 29, Issue 1, January 2005.

Use NaviTip tips to deliver these products and many others:

- **ChlorCid™ Surf** Sodium Hypochlorite Solution with Surfactant
- **ChlorCid™ & ChlorCid™ V** Sodium Hypochlorite Solutions (use with NaviTip™ 31 ga Double Sideport Irrigator Tip)
- **Ultradent™ EDTA 18% Solution**
- **Consepsis™ & Consepsis™ V** Chlorhexidine Antibacterial Solutions
- **File-Eze™** EDTA Lubricant
- **Ultradent™ Citric Acid 20% Solution**
- **UltraCal™ XS** Calcium Hydroxide Paste
- **EndoREZ™** Hydrophilic Resin Root Canal Sealer
- **MTAFlow™** Repair Cement

