

# NaviTip™ Delivery Tips

When it comes to root canals, Ultradent is with you every step of the way.

**NEW**



NaviTip™ Single Sideport Tip

NaviTip™ Double Sideport Tip

NaviTip™ Tip

NaviTip™ FX™ Tip

## Controlled Delivery

Ultradent's NaviTip tips provide the optimal irrigant delivery for every step in a root canal procedure. Unlike traditional irrigation tips, the cannula of each NaviTip tip is slightly rigid through the base and center, but flexible at the tip to allow for easy navigation of curved canals and precise delivery of irrigants short of the apex.

# NaviTip™



## NaviTip™ 29 ga or 30 ga Tip

- Controlled delivery anywhere in the canal
- Flexible cannula easily navigates curved canals

LOK-TITE™	Tip length	20pk	50pk
29 ga - 0.0129" NaviTip	27 mm	5115	1377
29 ga - 0.0129" NaviTip	25 mm	5114	1376
29 ga - 0.0129" NaviTip	21 mm	5113	1374
29 ga - 0.0129" NaviTip	17 mm	5112	1378
29 ga-29 ga NaviTips	27-17 mm	5116	1379
30 ga - 0.0118" NaviTip	27 mm	1354	1424
30 ga - 0.0118" NaviTip	25 mm	1250	1423
30 ga - 0.0118" NaviTip	21 mm	1349	1422
30 ga - 0.0118" NaviTip	17 mm	1249	1421
30 ga-30 ga NaviTips	27-17 mm	1351	3319



## NaviTip™ 29 ga Single Sideport Tip

- Designed to direct the flow of chemistry through the sideport of the tip before flowing down into the area of the apex, thus reducing the risk of product extrusion
- Flexible, stainless steel cannula easily navigates curved canals

LOK-TITE™	Tip length	20pk
29 ga - 0.0129" NaviTip	27 mm	4989
29 ga - 0.0129" NaviTip	25 mm	4990
29 ga - 0.0129" NaviTip	21 mm	4991
29 ga - 0.0129" NaviTip	17 mm	4992
29 ga-29 ga NaviTips	27-17 mm	5143

Note: UltraCal™ XS calcium hydroxide paste should only be used with NaviTip 29 ga Single Sideport tips.




## NaviTip™ 30 ga FX™ Tip

- Brush tip scrubs, cleans, and irrigates simultaneously
- Flexible cannula easily navigates curved canals
- Listed as an "EXCELLENT" product by a prominent independent research institute\*

\*Clinical Research Associates Newsletter, Volume 29, Issue 1, January 2005.

LOK-TITE™	Tip length	20pk
30 ga - 0.0118" NaviTip FX	17 mm	1452
30 ga - 0.0118" NaviTip FX	25 mm	1454

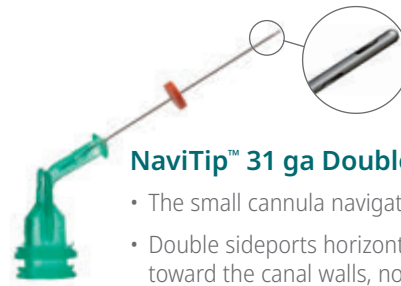

## NaviTip™ 31 ga Double Sideport Irrigator Tip


- The small cannula navigates the most intricate canals
- Double sideports horizontally express irrigant toward the canal walls, not the apex

LOK-TITE™	Tip length	20pk	50pk
31 ga - 0.0098" NaviTip	21 mm	5121	5122
31 ga - 0.0098" NaviTip	27 mm	5123	5124

## NaviTip™ Tips Adapt to Your Style

NaviTip tips deliver low-viscosity materials exactly where they are needed—even in the smallest canals.

The double sideport openings in the NaviTip 31 ga tip allow your irrigant to be safely expressed against the canal walls rather than toward—or even through—the foramen.

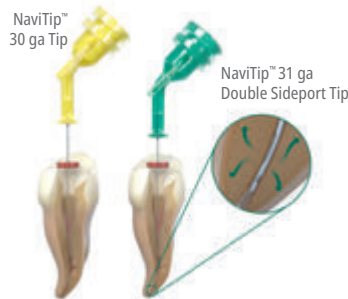

## Use NaviTip tips to deliver these products and many others:

- **ChlorCid™ Surf** Sodium Hypochlorite Solution with Surfactant
- **ChlorCid™ & ChlorCid™ V** Sodium Hypochlorite Solutions (use with NaviTip™ 31 ga Double Sideport Irrigator Tip)
- **Ultradent™ EDTA** 18% Solution
- **Consepsis™ & Consepsis™ V** Chlorhexidine Antibacterial Solutions

- **File-Eze™** EDTA Lubricant
- **Ultradent™ Citric Acid** 20% Solution
- **UltraCal™ XS** Calcium Hydroxide Paste (UltraCal™ XS calcium hydroxide paste should only be used with NaviTip 29 ga Single Sideport tips)
- **EndoREZ™** Hydrophilic Resin Root Canal Sealer
- **MTAFlow™** Repair Cements
- **MTApex™** Bioceramic Root Canal Sealer